

HPLC Application

ID No.: 10763

Dehydroacetic Acid in tap Water

Column: PhenoSphere-NEXT[™] 5 μ m C18 120 Å, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4308-E0

Elution Type: Isocratic

Eluent A: Acetonitrile

Eluent B: 0.1% Acetic acid

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	75	25

Flow Rate: 1 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 200 nm (ambient)

Analyst Note: SecurityGuard[™] Guard Cartridge System extends column lifetime.


- SecurityGuard Cartridges, C18 4 x 3.0mm, 10/Pk Part No.: AJ0-4287

- Holder Part No.: KJ0-4282



Products used in this application:

PhenoSphere[™]

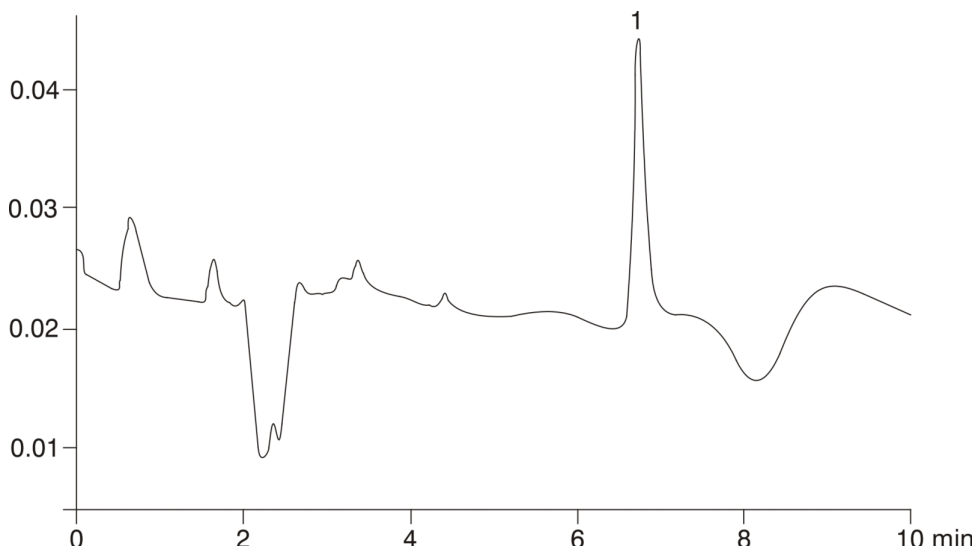

ANALYTES:

1 Dehydroacetic acid

