

# HPLC Application

ID No.: 1286

## Flavors

**Column:** Luna® 3 µm Phenyl-Hexyl 100 Å, LC Column 150 x 2 mm, Ea

**Dimensions:** 150 x 2 mm ID

**Order No:** 00F-4256-B0

**Elution Type:** Isocratic

**Eluent A:** 0.15% TFA in Water

**Eluent B:** Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	63	37

**Flow Rate:** 0.7 mL/min

**Col. Temp.:** 60 °C

**Detection:** Mass Spectrometer @ (ambient)


**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Phenyl 4 x 2.0mm, 10/Pk Part No.: AJ0-4350

- Holder Part No.: KJ0-4282



Products used in this application:

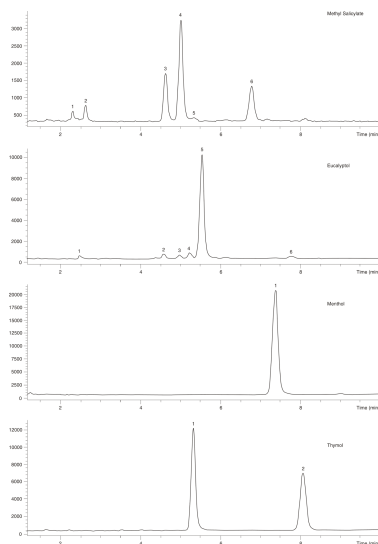

## ANALYTES:

- 1 Methyl salicylate
- 2 Eucalyptol
- 3 Menthol
- 4 Thymol

