

# HPLC Application


ID No.: 13604

## Chloroquine

**Column:** Chirex<sup>®</sup> (S)-LEU and (R)-NEA, LC Column 250 x 4.6 mm, Ea  
**Dimensions:** 250 x 4.6 mm ID  
**Order No:** 00G-3020-E0  
**Elution Type:** Isocratic  
**Eluent A:** Hexane/Dichloroethane/Ethanol-TFA (50:35:15) w/ ethanol-TFA (20:1)  
**Gradient**


Step No.	Time (min)	Pct A
1	0	100

**Profile:**  
**Flow Rate:** 1 mL/min  
**Col. Temp.:** ambient  
**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 346 nm (ambient)

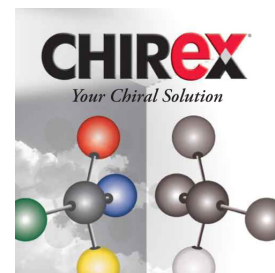

Products used in this application:

**CHIREX**<sup>®</sup>

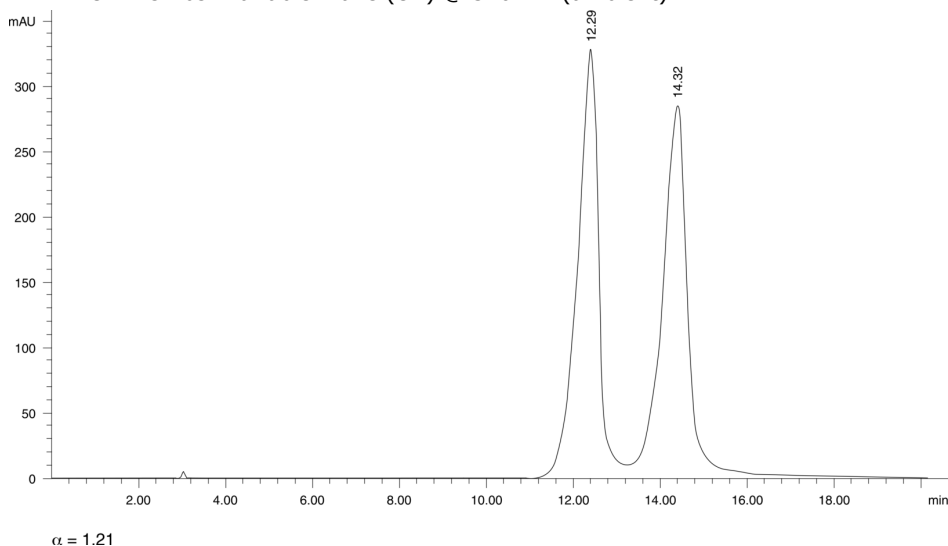

### ANALYTES:

- 1 Chloroquine  
Retention Time: 12.29 min
- 2 Chloroquine  
Retention Time: 14.32 min

