

# HPLC Application


ID No.: 13628

## Doxylamine

**Column:** Chirex® (S)-VAL and (R)-NEA, LC Column 250 x 4.6 mm, Ea  
**Dimensions:** 250 x 4.6 mm ID  
**Order No:** 00G-3014-E0  
**Elution Type:** Isocratic  
**Eluent A:** Hexane/Dichloroethane/Ethanol-TFA 55:35:10 ethanol-TFA 20:1  
**Gradient**


Step No.	Time (min)	Pct A
1	0	100

**Profile:**  
**Flow Rate:** 0.7 mL/min  
**Col. Temp.:** ambient  
**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 262 nm (ambient)

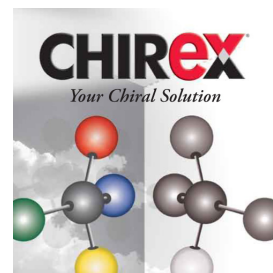

Products used in this application:

**CHIREX**®

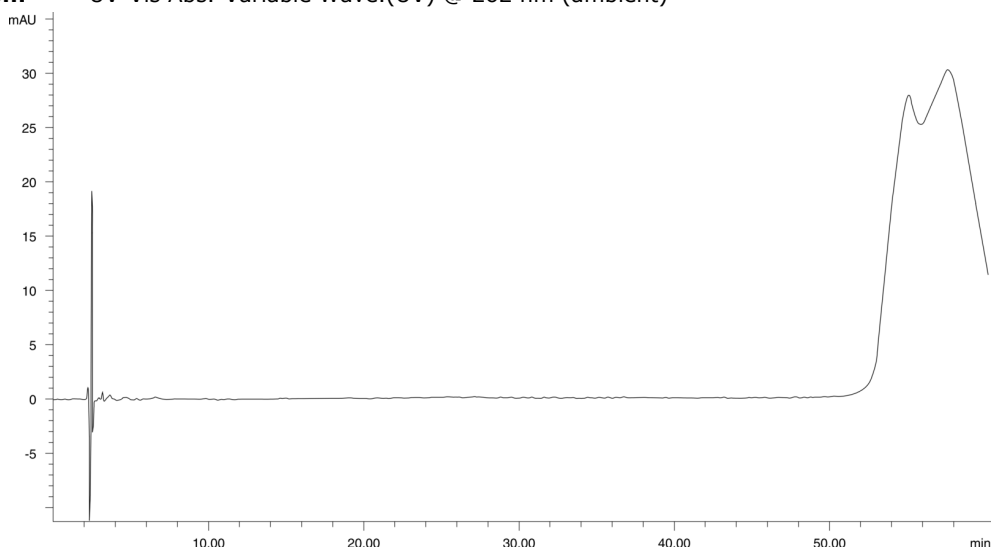

### ANALYTES:

1 Doxylamine

