

# HPLC Application


ID No.: 13659

## Methoxamine

**Column:** Chirex® (S)-VAL and (R)-NEA, LC Column 250 x 4.6 mm, Ea  
**Dimensions:** 250 x 4.6 mm ID  
**Order No:** 00G-3014-E0  
**Elution Type:** Isocratic  
**Eluent A:** Hexane/Dichloroethane/Ethanol-TFA (76:20:4) ethanol-TFA 20:1  
**Gradient**


Step No.	Time (min)	Pct A
1	0	100

**Profile:**  
**Flow Rate:** 0.7 mL/min  
**Col. Temp.:** ambient  
**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 290 nm (ambient)

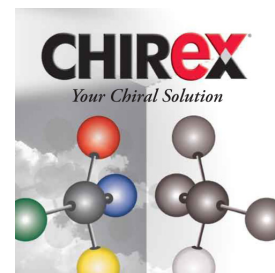

Products used in this application:

**CHIREX**®

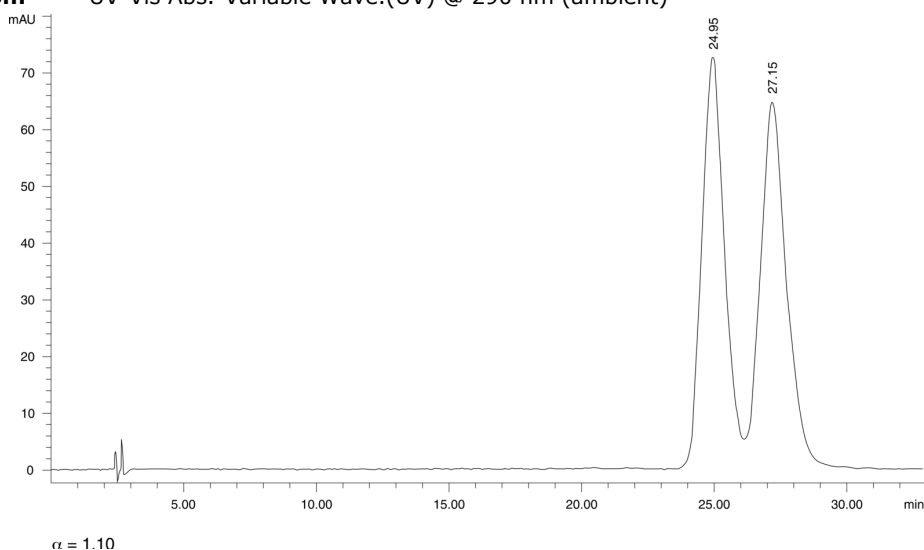

### ANALYTES:

- 1 Methoxamine  
Retention Time: 24.95 min
- 2 Methoxamine  
Retention Time: 27.15 min

