

# HPLC Application

ID No.: 14041

## Ornithine

**Column:** Chirex® 3126 (D)-penicillamine, LC Column 150 x 4.6 mm, Ea

**Dimensions:** 150 x 4.6 mm ID

**Order No:** 00F-3126-E0

**Elution Type:** Isocratic

**Eluent A:** 1 mM Copper (II) sulfate in water

**Gradient**



Step No.	Time (min)	Pct A
1	0	100

**Profile:**

**Flow Rate:** 0.2 mL/min

**Col. Temp.:** ambient

**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)

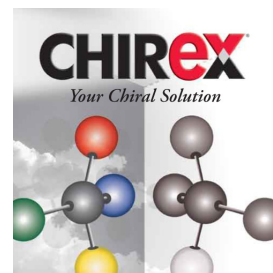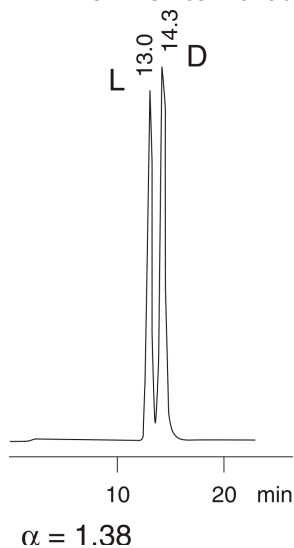

Products used in this application:

**CHIREX**®

### ANALYTES:

1 L-Ornithine

Retention Time: 13.833 min

2 D-Ornithine

Retention Time: 14.283 min

