

# HPLC Application

ID No.: 14089

## 1-Aminoethylphosphonic acid

**Column:** Chirex<sup>®</sup> 3126 (D)-penicillamine, LC Column 150 x 4.6 mm, Ea

**Dimensions:** 150 x 4.6 mm ID

**Order No:** 00F-3126-E0

**Elution Type:** Isocratic

**Eluent A:** 2 mM Copper (II) sulfate in water


**Eluent B:** Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	95	5

**Flow Rate:** 1 mL/min

**Col. Temp.:** ambient

**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)

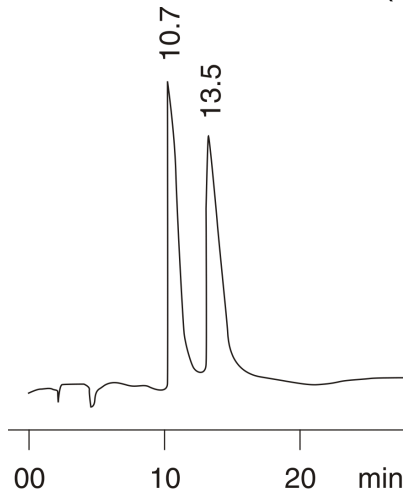

$\alpha = 1.30$


### ANALYTES:

**1** 1-Aminoethylphosphonic acid

Retention Time: 10.733 min

**2** 1-Aminoethylphosphonic acid

Retention Time: 13.548 min

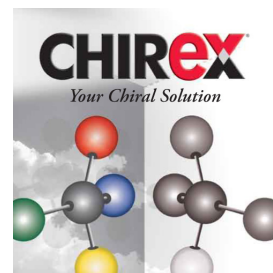

Products used in this application:

**CHIREX**<sup>®</sup>

