

# HPLC Application

ID No.: 14113

## 3,4-Dihydro-2H-pyran-2-carboxylic acid

**Column:** Chirex<sup>®</sup> 3126 (D)-penicillamine, LC Column 150 x 4.6 mm, Ea

**Dimensions:** 150 x 4.6 mm ID

**Order No:** 00F-3126-E0

**Elution Type:** Isocratic

**Eluent A:** 2 mM Copper (II) sulfate in water

**Eluent B:** Isopropanol

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	90	10

**Flow Rate:** 1 mL/min


**Col. Temp.:** ambient

**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)



Products used in this application:

**CHIREX**<sup>®</sup>

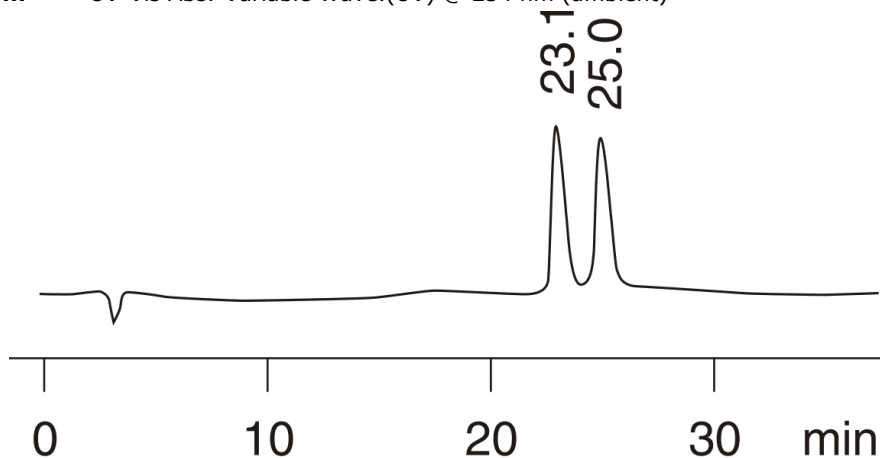

$$\alpha = 1.09$$

### ANALYTES:

- 1 3,4-Dihydro-2H-pyran-2-carboxylic acid

Retention Time: 23.067 min

- 2 3,4-Dihydro-2H-pyran-2-carboxylic acid

Retention Time: 25.033 min

