

# HPLC Application

ID No.: 14253

## Opiates on Synergi Max-RP

**Column:** Synergi<sup>™</sup> 4 µm Max-RP 80 Å, LC Column 150 x 4.6 mm, Ea

**Dimensions:** 150 x 4.6 mm ID

**Order No:** 00F-4337-E0

**Elution Type:** Gradient

**Eluent A:** 20mM Potassium dihydrogen phosphate

**Eluent B:** Acetonitrile

<b>Gradient Profile:</b>	<b>Step No.</b>	<b>Time (min)</b>	<b>Pct A</b>	<b>Pct B</b>
	<b>1</b>	0	95	5
	<b>2</b>	15	65	35

**Flow Rate:** 1.5 mL/min


**Col. Temp.:** ambient

**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 210 nm (ambient)

**Analyst Note:** SecurityGuard<sup>™</sup> Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, C12 4 x 3.0mm, 10/Pk Part No.: AJ0-6074

- Holder Part No.: KJ0-4282

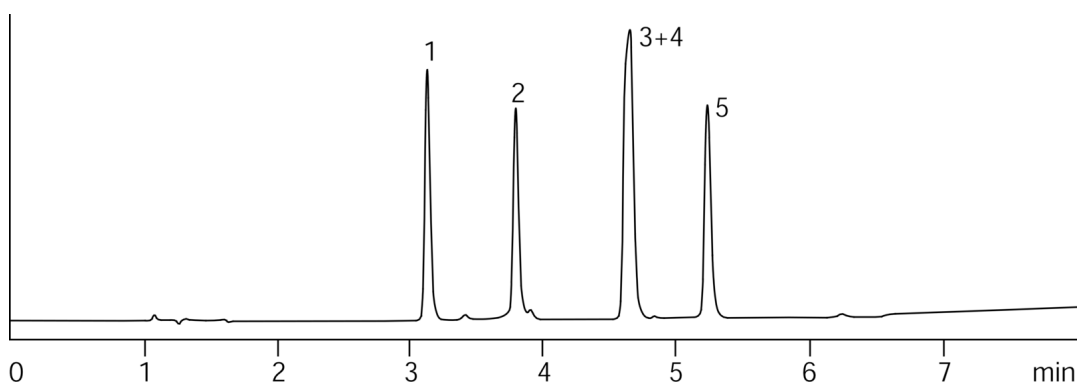

### ANALYTES:

- 1 Morphine
- 2 Hydrpmorphone
- 3 Codeine / Nalorphine
- 4 Oxycodone

