

## Acidic compounds on Luna Phenyl-Hexyl

**Column:** Luna® 5 µm Phenyl-Hexyl 100 Å, LC Column 150 x 4.6 mm, Ea

**Dimensions:** 150 x 4.6 mm ID

**Order No:** 00F-4257-E0

**Elution Type:** Isocratic

**Eluent A:** 20mM Potassium phosphate, pH 2.5

**Eluent B:** Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	75	25

**Flow Rate:** 1.5 mL/min

**Col. Temp.:** 30 °C

**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)


**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Phenyl 4 x 3.0mm, 10/Pk Part No.: AJ0-4351

- Holder Part No.: KJ0-4282



Products used in this application:

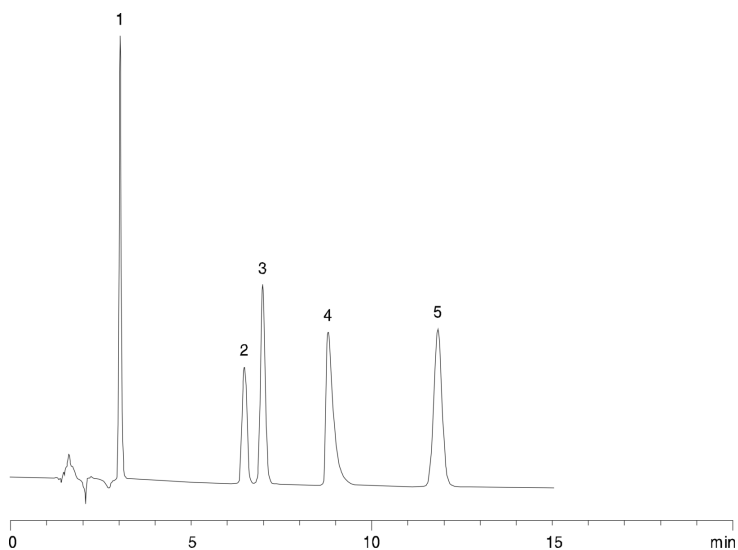

### ANALYTES:

- 1 p-Hydroxybenzoic acid
- 2 Sorbic acid
- 3 Benzoic acid
- 4 Salicylic acid
- 5 p-Toluic acid

