

HPLC Application

ID No.: 14467

EPA Method 610 - PAH's on Luna 5um C18(2)

Column: Luna® 5 µm C18(2) 100 Å, LC Column 250 x 4.6 mm, Ea

Dimensions: 250 x 4.6 mm ID

Order No: 00G-4252-E0

Elution Type: Gradient

Eluent A: Water

Eluent B: Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	60	40
	2	5	60	40
	3	30	0	100

Flow Rate: 1.55 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 210 nm (ambient)


Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, C18 4 x 3.0mm, 10/Pk Part No.: AJ0-4287

- Holder Part No.: KJ0-4282



Products used in this application:

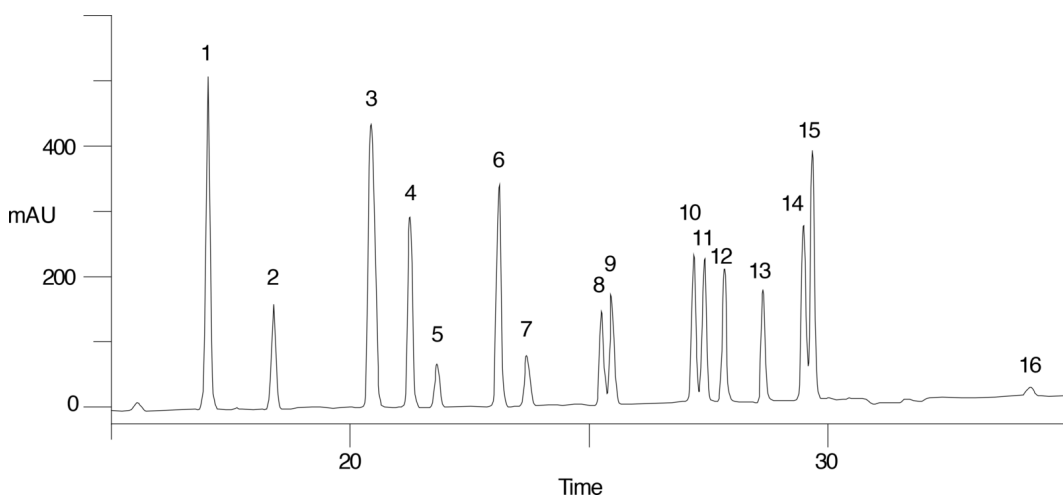

EPA Method 610 - PAH's on Luna 5um C18(2)

ANALYTES:

- 1 Naphthalene
- 2 Acenaphthylene
- 3 Acenaphthene
- 4 Fluorene
- 5 Phenanthrene
- 6 Anthracene
- 7 Fluoranthene
- 8 Pyrene
- 9 Benz[a]anthracene
- 10 Chrysene
- 11 Benzo[b]fluoranthene
- 12 Benzo[k]fluoranthene
- 13 Benzo[a]pyrene
- 14 Dibenzo[a,h]anthracene
- 15 Benzo[g,h,i]perylene
- 16 Indeno[1,2,3-cd]pyrene

