


## 12-18mer Poly-dT Oligonucleotide Mix On Jupiter Proteo

**Column:** N/A  
**Dimensions:** N/A  
**Order No:** N/A  
**Elution Type:** N/A  
**Eluent A:** 0.1M TEAA, pH 7.0  
**Eluent B:** Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	97	3
	2	20	85	15

**Flow Rate:** 1 mL/min  
**Col. Temp.:** 60 °C  
**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 260 nm (ambient)  
**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, C12 4 x 3.0mm, 10/Pk Part No.: AJ0-6074
- Holder Part No.: KJ0-4282

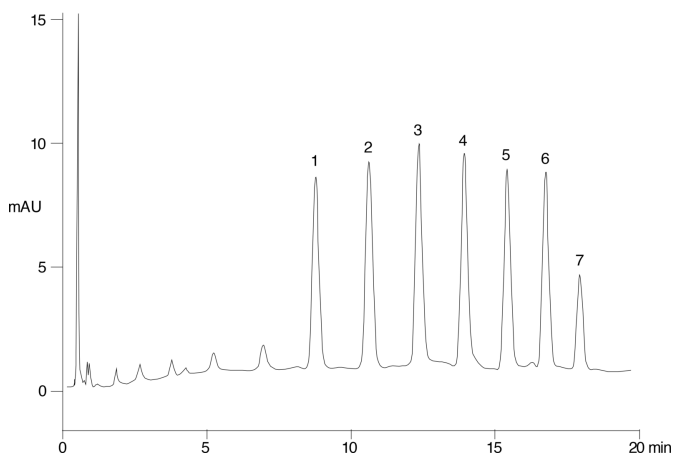

### ANALYTES:

- 1 12-18mer Poly-dT oligonucleotide mixture

