

12-18mer Poly-dT Oligonucleotide Mix On Jupiter Proteo

Column: N/A
Dimensions: N/A
Order No: N/A
Elution Type: N/A
Eluent A: 0.1M TEAA, pH 7.0
Eluent B: Acetonitrile

| Gradient Profile: | Step No. | Time (min) | Pct A | Pct B |
|-------------------|----------|------------|-------|-------|
| | 1 | 0 | 97 | 3 |
| | 2 | 20 | 85 | 15 |

Flow Rate: 1 mL/min


Col. Temp.: 60 °C

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 260 nm (ambient)

Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, C12 4 x 3.0mm, 10/Pk Part No.: AJ0-6074

- Holder Part No.: KJ0-4282

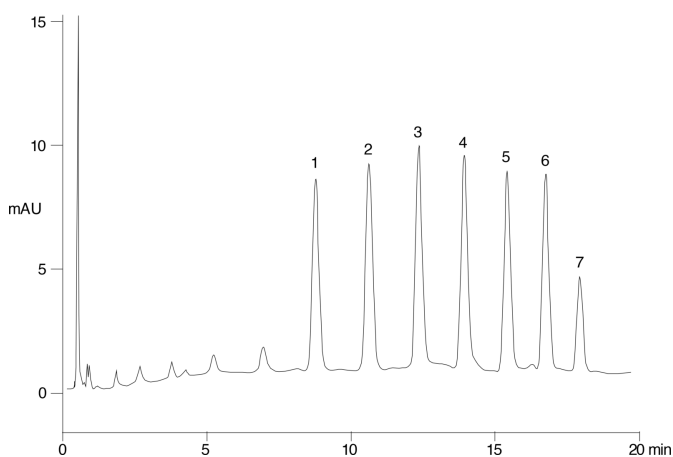

ANALYTES:

1 12-18mer Poly-dT oligonucleotide mixture

