


12-18mer Poly-dT Oligonucleotide Mix On Jupiter Proteo

Column: N/A
Dimensions: N/A
Order No: N/A
Elution Type: N/A
Eluent A: 0.1M TEAA, pH 7.0
Eluent B: Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	97	3
	2	20	85	15

Flow Rate: 1 mL/min
Col. Temp.: 60 °C
Detection: UV-Vis Abs.-Variable Wave.(UV) @ 260 nm (ambient)
Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, C12 4 x 3.0mm, 10/Pk Part No.: AJ0-6074
- Holder Part No.: KJ0-4282

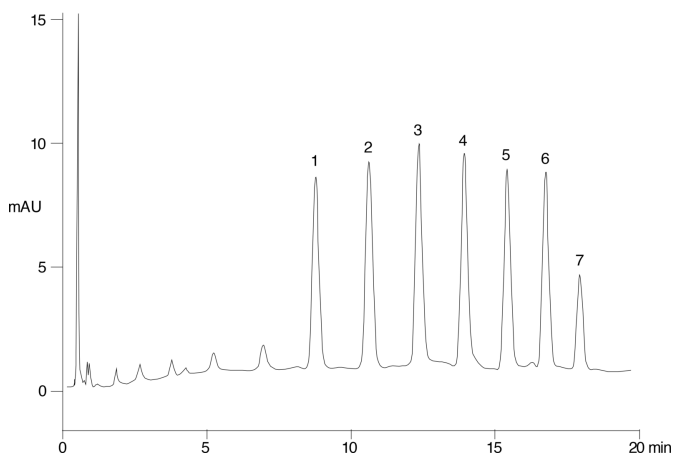

ANALYTES:

- 1 12-18mer Poly-dT oligonucleotide mixture

