

## Tricyclic Antidepressants on Hyperclone 5u C8 BDS

**Column:** HyperClone™ 5 µm BDS C8 130 Å, LC Column 150 x 4.6 mm, Ea

**Dimensions:** 150 x 4.6 mm ID

**Order No:** 00F-4418-E0

**Elution Type:** Isocratic

**Eluent A:** 20mM KH<sub>2</sub>PO<sub>4</sub>, pH=7.0

**Eluent B:** CH<sub>3</sub>CN

**Eluent C:** MeOH

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B	Pct C
	1	0	30	35	35

**Flow Rate:** 1 mL/min

**Col. Temp.:** ambient

**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)

**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.


- SecurityGuard Cartridges, C8 4 x 3.0mm, 10/Pk Part No.: AJ0-4290

- Holder Part No.: KJ0-4282



Products used in this application:

**HyperClone™**

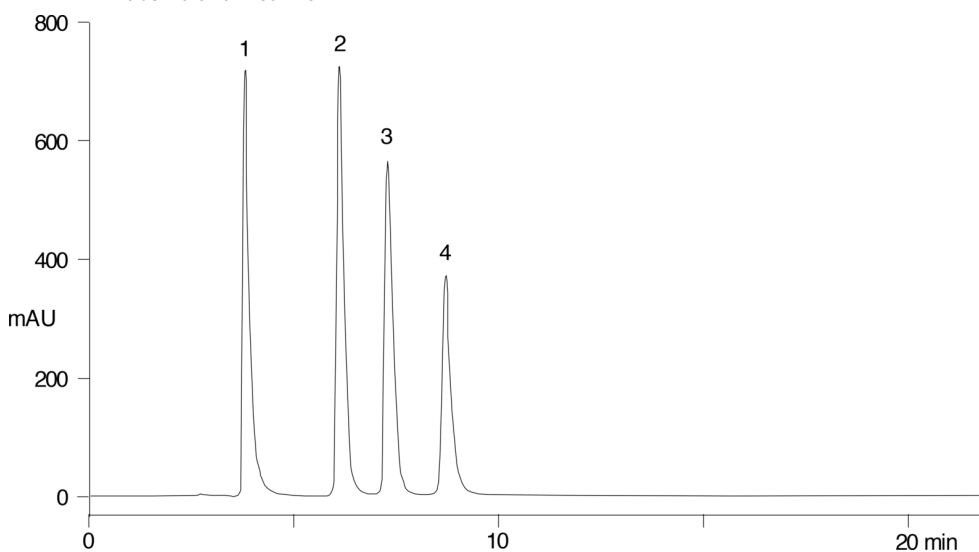

### ANALYTES:

- 1 Imipramine
- 2 Nortriptyline
- 3 Amitriptyline
- 4 Clomipramine

