

Industrial Solvents on ZB-624 60x0.32x1.8

Column: Zebtron™ ZB-624, GC Cap. Column 60 m x 0.32 mm x 1.80 μm, Ea
Phase: 6% Cyanopropylphenyl 94% Dimethylpolysiloxane
Dimensions: 60 meters x 0.32 mm x 1.8 μm
Order No: 7KM-G005-31
Oven Profile: 40°C for 8 min. to 90°C at 8°C/min. to 210°C at 15°C for 8 min

Carrier Gas: Constant Flow Helium, 1 mL/min

Injection: Split 20:1 0.6 μL @ 220°C

Detection: Mass Selective (MSD) (180°C)


Analyst Note: All analytes were 140ng on column except methylene chloride as the solvent.



Products used in this application:



14807

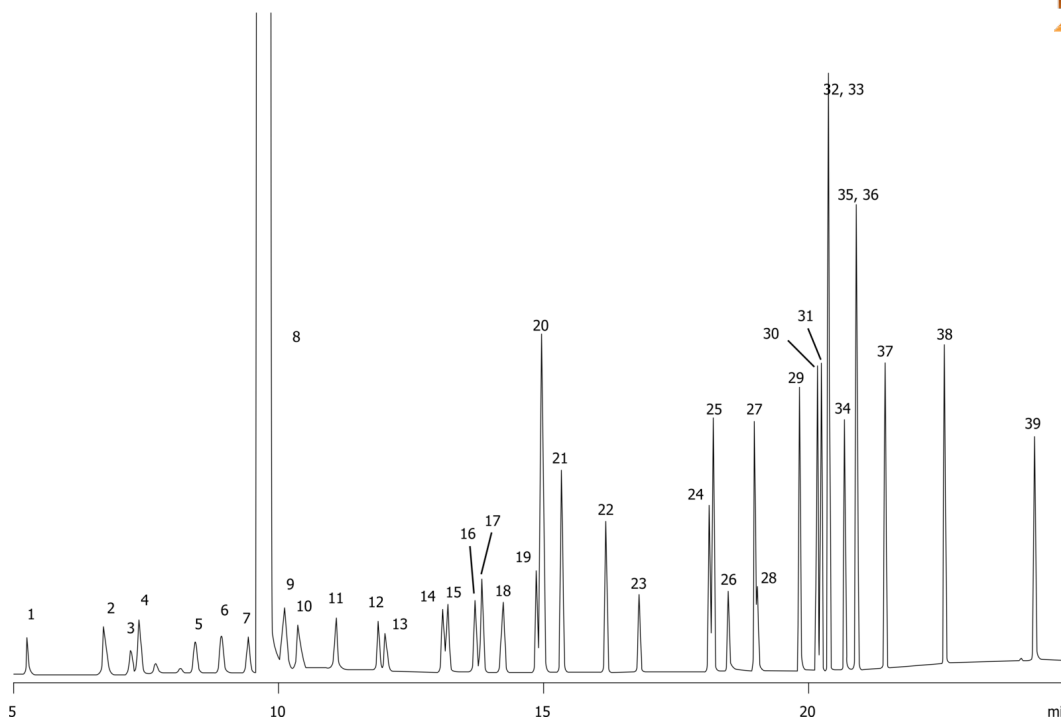

Industrial Solvents on ZB-624 60x0.32x1.8

ANALYTES:

- 1 Methanol
- 2 Pentane
- 3 Ethanol
- 4 Ether
- 5 Acetone
- 6 Isopropanol
- 7 Acetonitrile
- 8 Methylene chloride
- 9 tert-Butanol
- 10 Methyl tert-butyl ether
- 11 Hexane
- 12 Allyl alcohol
- 13 n-Propanol
- 14 Methyl ethyl ketone
- 15 Ethyl acetate
- 16 Tetrahydrofuran
- 17 Chloroform
- 18 Cyclohexane
- 19 1,2-Dimethoxyethane
- 20 Benzene
- 21 Heptane
- 22 Trichloroethylene
- 23 1,4-Dioxane
- 24 Pyridine
- 25 Toluene
- 26 Ethylene glycol
- 27 Tetrachloroethylene
- 28 Formamide
- 29 N,N-Dimethylformamide
- 30 Chlorobenzene
- 31 Ethylbenzene
- 32 m-Xylene
- 33 p-Xylene
- 34 3-Picoline
- 35 o-Xylene
- 36 Styrene
- 37 Anisole
- 38 Benzaldehyde
- 39 N-Methylpyrrolidone

