

# HPLC Application

ID No.: 14847

## Basic drugs at pH 9 using Synergi Fusion-RP

**Column:** Synergi<sup>™</sup> 4 µm Fusion-RP 80 Å, LC Column 150 x 4.6 mm, Ea

**Dimensions:** 150 x 4.6 mm ID

**Order No:** 00F-4424-E0

**Elution Type:** Gradient

**Eluent A:** 25mM TRIS, pH 9

**Eluent B:** Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0.5	60	40
	2	15	30	70

**Flow Rate:** 1 mL/min

**Col. Temp.:** 30 °C

**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)


**Analyst Note:** SecurityGuard<sup>™</sup> Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Fusion-RP 4 x 3.0mm, 10/Pk Part No.: AJ0-7557

- Holder Part No.: KJ0-4282



14847

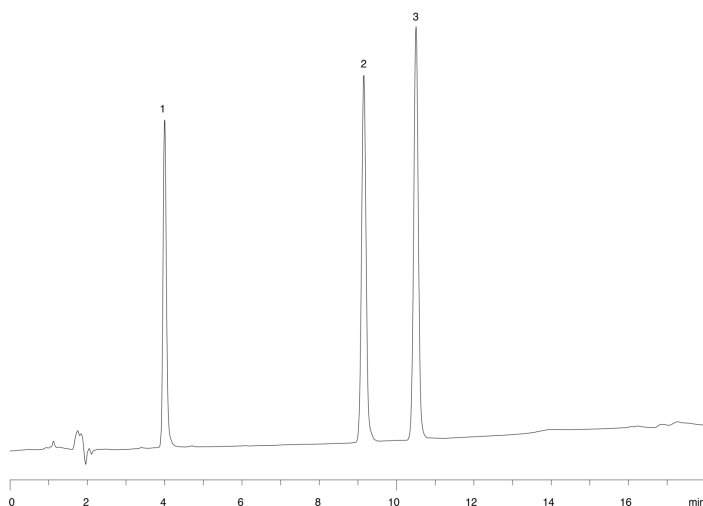

### ANALYTES:

- 1 Scopolamine
- 2 Lidocaine
- 3 Noscapine

