


Nitrosamines ZB-WAX 30x0.32x0.5

Column: Zebtron™ ZB-WAX, GC Cap. Column 30 m x 0.32 mm x 0.50 µm, Ea
Phase: Polyethylene Glycol (PEG)
Dimensions: 30 meters x 0.32 mm x 0.5 µm
Order No: 7HM-G007-17
Oven Profile: 40°C to 260°C at 20°C/min for 5 min.

Carrier Gas: Constant Flow Helium, 1.4 mL/min
Injection: split 50:1 1 µL @ 200°C
Detection: Mass Selective (MSD) (240°C)

14854

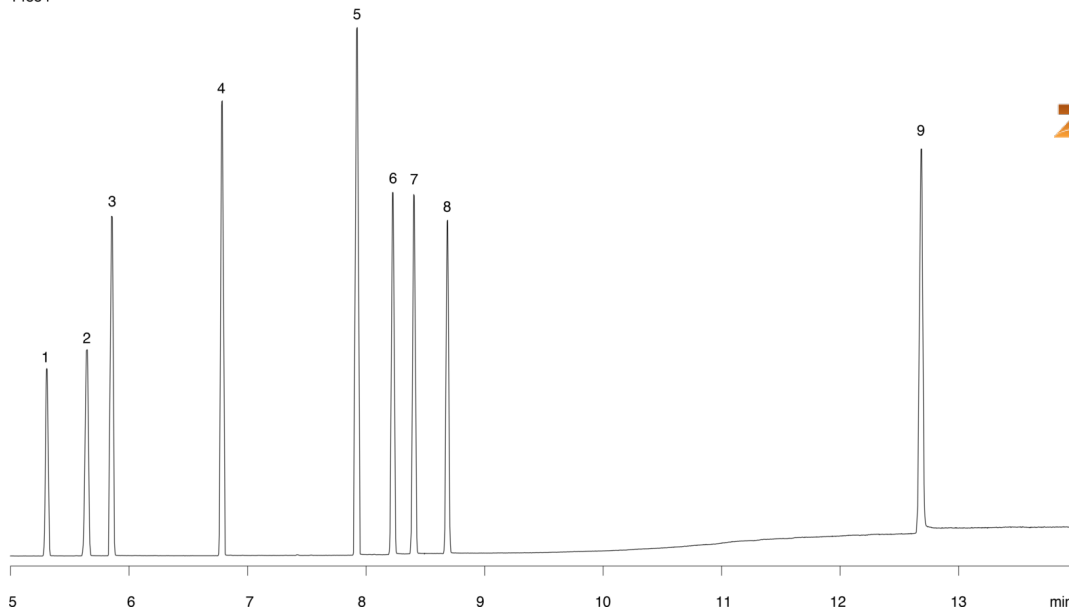

Products used in this application:



ANALYTES:

- 1 N-Nitrosodimethylamine
- 2 N-Nitrosomethylethylamine
- 3 N-Nitrosodiethylamine
- 4 N-Nitrosodi-n-propylamine
- 5 N-Nitrosodi-n-butylamine
- 6 N-Nitrosopiperidine
- 7 N-Nitrosopyrrolidine
- 8 N-Nitrosomorpholine
- 9 N-Nitrosodiphenylamine

