

## Nitrosamines ZB-WAX 30x0.32x0.5

**Column:** Zebtron™ ZB-WAX, GC Cap. Column 30 m x 0.32 mm x 0.50 µm, Ea

**Phase:** Polyethylene Glycol (PEG)

**Dimensions:** 30 meters x 0.32 mm x 0.5 µm

**Order No:** 7HM-G007-17


**Oven Profile:** 40°C to 260°C at 20°C/min for 5 min.

**Carrier Gas:** Constant Flow Helium, 1.4 mL/min

**Injection:** split 50:1 1 µL @ 200°C

**Detection:** Mass Selective (MSD) (240°C)

14854

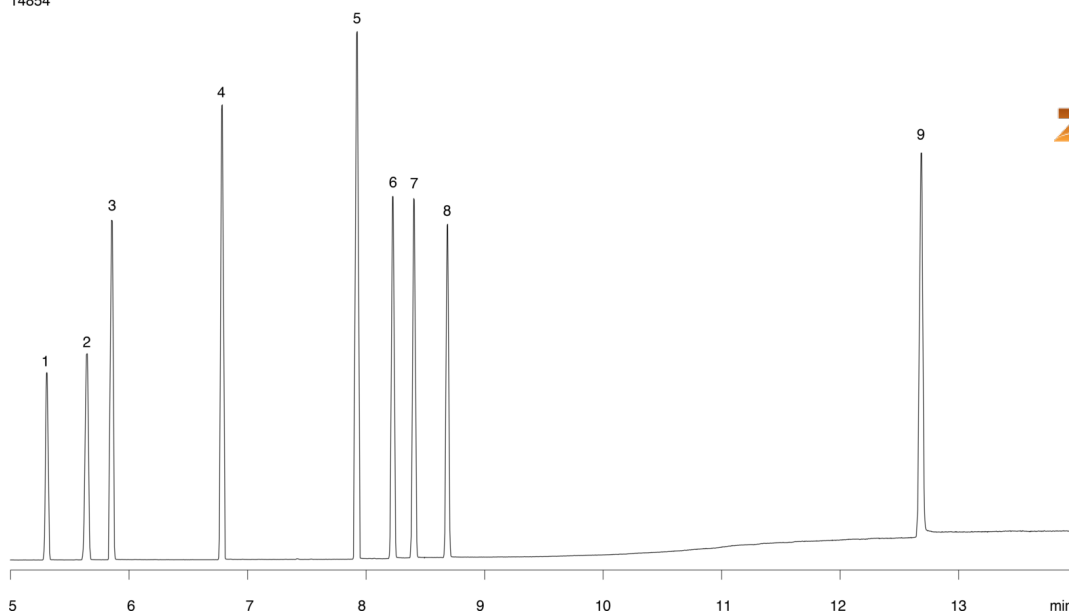

Products used in this application:



### ANALYTES:

- 1 N-Nitrosodimethylamine
- 2 N-Nitrosomethylethylamine
- 3 N-Nitrosodiethylamine
- 4 N-Nitrosodi-n-propylamine
- 5 N-Nitrosodi-n-butylamine
- 6 N-Nitrosopiperidine
- 7 N-Nitrosopyrrolidine
- 8 N-Nitrosomorpholine
- 9 N-Nitrosodiphenylamine

