

# HPLC Application

ID No.: 15455

## Benzodiazepines (SPE-Plasma Blank) on Gemini C18

**Column:** Gemini<sup>®</sup> 5  $\mu$ m C18 110 Å, LC Column 150 x 3 mm, Ea

**Dimensions:** 150 x 3 mm ID

**Order No:** 00F-4435-Y0

**Elution Type:** Gradient

**Eluent A:** 0.1% Formic acid in Water

**Eluent B:** 0.1% Formic acid in Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	80	20
	2	8	20	80
	3	9	20	80

**Flow Rate:** 0.8 mL/min

**Col. Temp.:** ambient

**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)

**Analyst Note:** SecurityGuard<sup>™</sup> Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Gemini C18 4 x 2.0mm, 10/Pk Part No.: AJ0-7596


- Holder Part No.: KJ0-4282



Products used in this application:



15455

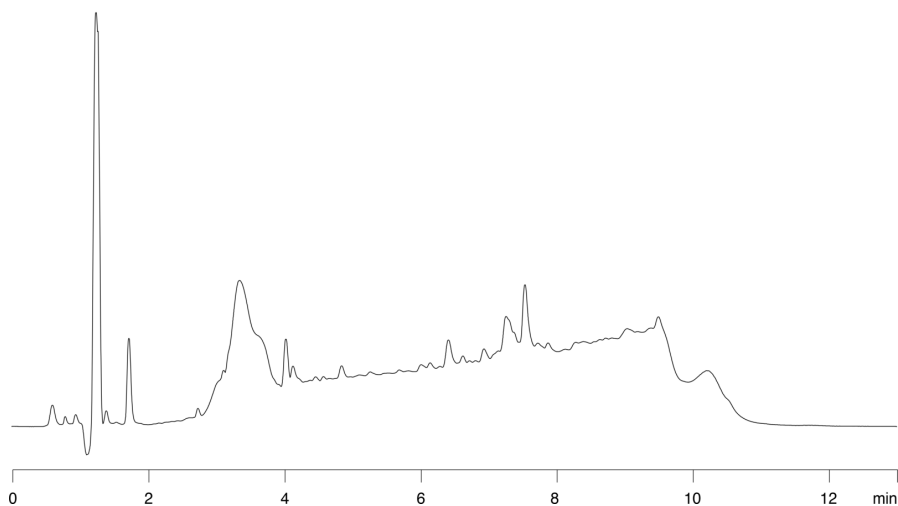

### ANALYTES:

1 blank

