

# HPLC Application

ID No.: 15460

## Selectivity of Acidic, Neutral and Basic Analytes at pH10.0

**Column:** Gemini<sup>®</sup> 5  $\mu$ m C18 110 Å, LC Column 150 x 4.6 mm, Ea

**Dimensions:** 150 x 4.6 mm ID

**Order No:** 00F-4435-E0

**Elution Type:** Isocratic

**Eluent A:** 10mM Ammonium Bicarbonate, pH 10

**Eluent B:** ACN

| Gradient Profile: | Step No. | Time (min) | Pct A | Pct B |
|-------------------|----------|------------|-------|-------|
|                   | 1        | 0          | 50    | 50    |

**Flow Rate:** 1 mL/min

**Col. Temp.:** ambient

**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)

**Analyst Note:** SecurityGuard<sup>™</sup> Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Gemini C18 4 x 3.0mm, 10/Pk Part No.: AJ0-7597


- Holder Part No.: KJ0-4282



Products used in this application:



15460

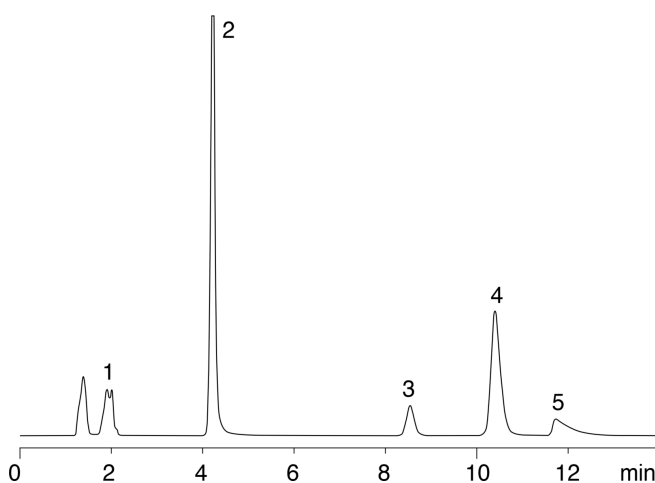

### ANALYTES:

- 1 Chlorpropamide
- 2 Butyl paraben
- 3 Lidocaine
- 4 Triprolidine
- 5 Dextromethorphan

