


## Extended List Polynuclear Aromatic Hydrocarbons (PAHs) on ZB-5MS

**Column:** Zebtron™ ZB-5MS, GC Cap. Column 20 m x 0.18 mm x 0.18 μm, Ea  
**Phase:** 5% Phenyl-Arylene 95% Dimethylpolysiloxane  
**Dimensions:** 20 meters x 0.18 mm x 0.18 μm  
**Order No:** 7FD-G010-08  
**Oven Profile:** 130°C for 0.5min to 250°C at 25°C/min to 270°C at 6°C/min to 320°C at 25°C/min for 4 min  
**Carrier Gas:** Constant Flow Helium, 1.2 mL/min  
**Injection:** Split 15:1 1 μL @ 285°C  
**Detection:** Mass Selective (MSD) (180°C)

15562



Products used in this application:

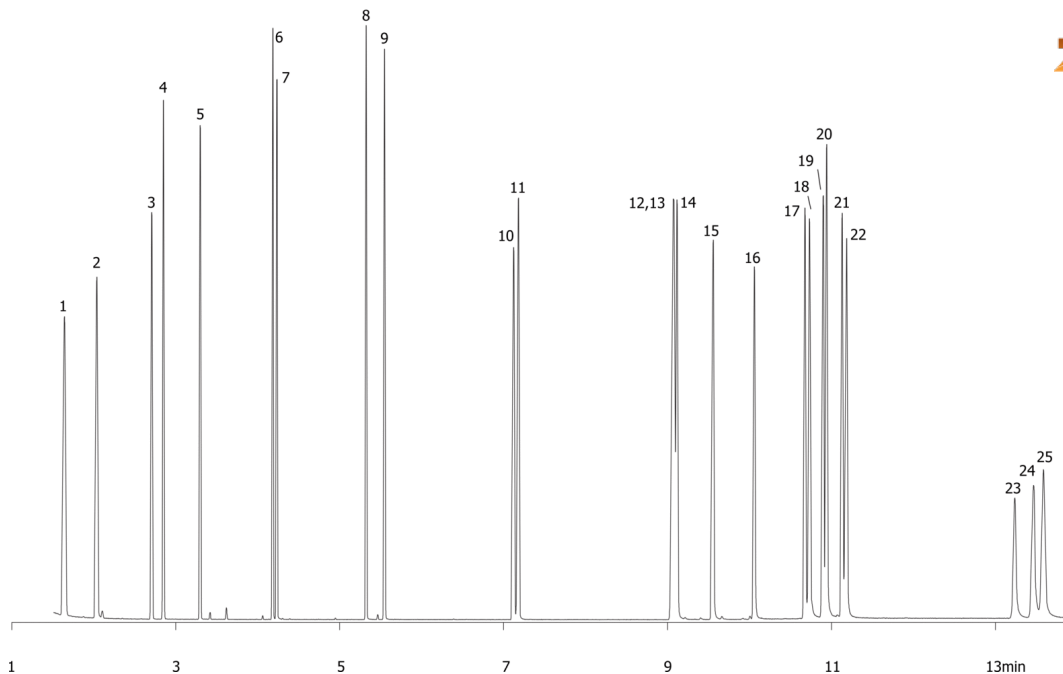

## Extended List Polynuclear Aromatic Hydrocarbons (PAHs) on ZB-5MS

### ANALYTES:

- 1 Naphthalene
- 2 2-Methylnaphthalene
- 3 Acenaphthalene
- 4 Acenaphthene
- 5 Fluorene
- 6 Phenanthrene
- 7 Anthracene
- 8 Fluoranthene
- 9 Pyrene
- 10 Benz[a]anthracene
- 11 Chrysene
- 12 Benzo[b]fluoranthene
- 13 Benzo[j]fluoranthene
- 14 Benzo[k]fluoranthene
- 15 Benzo[a]pyrene
- 16 3-Methylcholanthrene
- 17 Dibenz[a,h]acridine
- 18 Dibenz[a,,j]acridine
- 19 Indeno[1,2,3-cd]pyrene
- 20 Dibenz[a,h]anthracene
- 21 7H-Dibenzo[c,g]carbazole
- 22 Benzo[g,h,i]perylene
- 23 Dibenzo[a,e]pyrene
- 24 Dibenzo[a,i]pyrene
- 25 Dibenzo[a,j]pyrene

