

# HPLC Application

ID No.: 15666

## Hydrocortisone USP modified (Synergi 4µm Fusion-RP)

**Column:** Synergi<sup>™</sup> 4 µm Fusion-RP 80 Å, LC Column 150 x 4.6 mm, Ea

**Dimensions:** 150 x 4.6 mm ID

**Order No:** 00F-4424-E0

**Elution Type:** Isocratic

**Eluent A:** Water

**Eluent B:** Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	75	25

**Flow Rate:** 1.3 mL/min


**Col. Temp.:** 25 °C

**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (25 °C)

**Analyst Note:** SecurityGuard<sup>™</sup> Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Fusion-RP 4 x 3.0mm, 10/Pk Part No.: AJ0-7557

- Holder Part No.: KJ0-4282

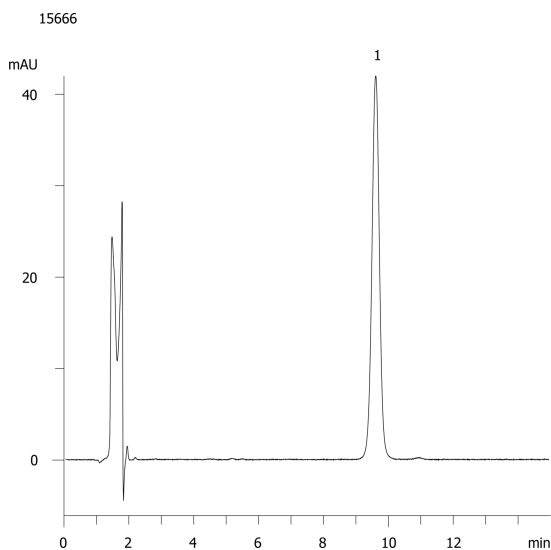

### ANALYTES:

1 Hydrocortisone

Retention Time: 9.61 min

