

Benzodiazepines on ZB-1MS

Column: Zebtron[™] ZB-1MS, GC Cap. Column 12 m x 0.20 mm x 0.33 µm, Ea

Phase: 100% Dimethylpolysiloxane

Dimensions: 12 meters x 0.2 mm x 0.33 µm

Order No: 7DE-G011-14

Oven Profile: 150°C to 200°C at 30°C/min to 300 at 15°C/min for 2 min

Carrier Gas: Constant Flow Helium, 1.2 mL/min


Injection: Splitless :1 3 µL @ 280°C

Detection: Mass Selective (MSD) (0°C)

15782



Products used in this application:

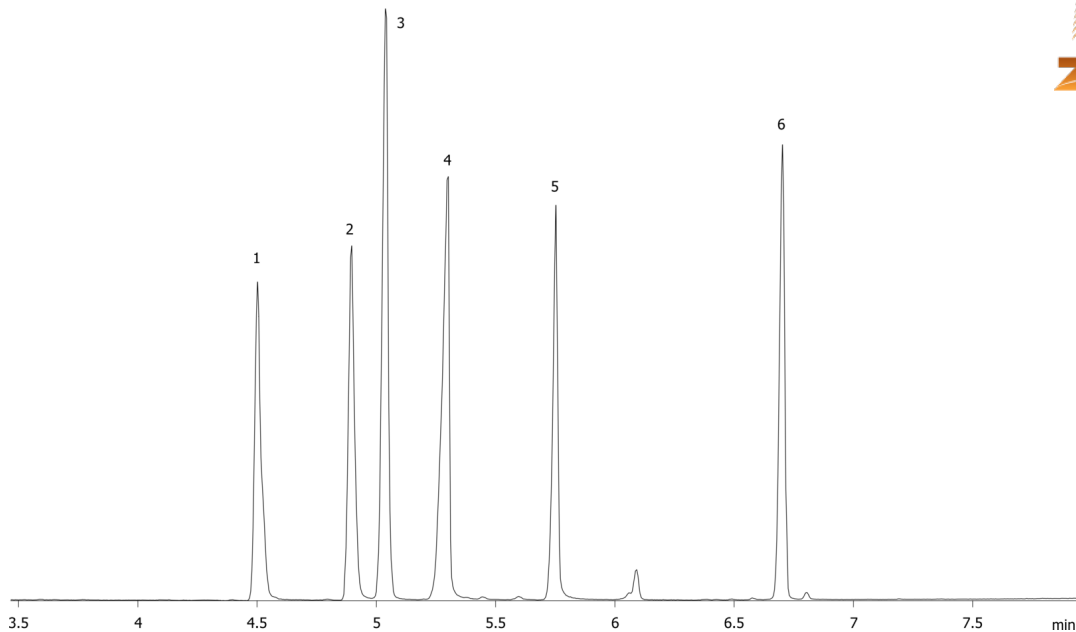

ANALYTES:

- 1 Oxazepam
- 2 Lorazepam
- 3 Diazepam
- 4 Nordazepam
- 5 Temazepam
- 6 Flurazepam

