

# HPLC Application

ID No.: 15858

## Non-Polar Basics at Low pH on Gemini 5um C6-Phenyl

**Column:** Gemini<sup>®</sup> 5  $\mu$ m C6-Phenyl 110  $\text{\AA}$ , LC Column 150 x 4.6 mm, Ea

**Dimensions:** 150 x 4.6 mm ID

**Order No:** 00F-4444-E0

**Elution Type:** Isocratic

**Eluent A:** 20mM Phosphate buffer pH2.5

**Eluent B:** Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	75	25

**Flow Rate:** 1.5 mL/min

**Col. Temp.:** ambient

**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 210 nm (ambient)


**Analyst Note:** SecurityGuard<sup>™</sup> Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Gemini C6-Phenyl 4 x 3.0mm, 10/Pk Part No.: AJ0-7915

- Holder Part No.: KJ0-4282



Products used in this application:

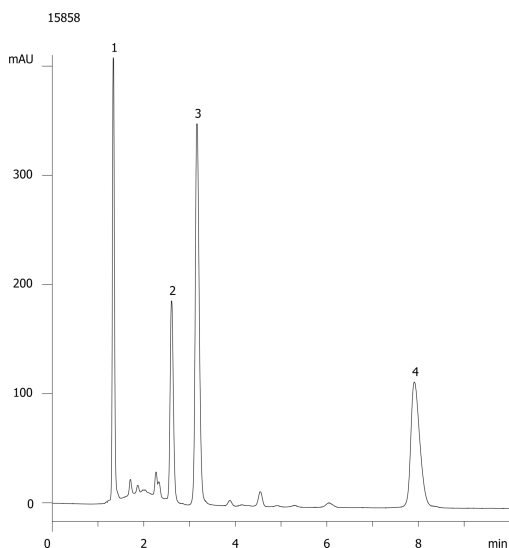

### ANALYTES:

- 1 Maleic acid
- 2 Chlorpheniramine
- 3 Triprolidine
- 4 Diphenhydramine

