

# HPLC Application

ID No.: 15949

## 40nt DNA Desalting with 25mM TBAC - Load, Wash, Elution Fractions

**Column:** Clarity<sup>®</sup> 3  $\mu$ m Oligo-RP, LC Column 50 x 4.6 mm, Ea

**Dimensions:** 50 x 4.6 mm ID

**Order No:** 00B-4441-E0

**Elution Type:** Gradient

**Eluent A:** 50mM TEAA/5% ACN

**Eluent B:** Methanol

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	90	10
	2	20	40	60

**Flow Rate:** 1 mL/min

**Col. Temp.:** ambient

**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 260 nm (ambient)


**Analyst Note:** SecurityGuard<sup>™</sup> Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Clarity Oligo-RP 4 x 3.0mm, 10/Pk Part No.: AJ0-8135

- Holder Part No.: KJ0-4282



Products used in this application:

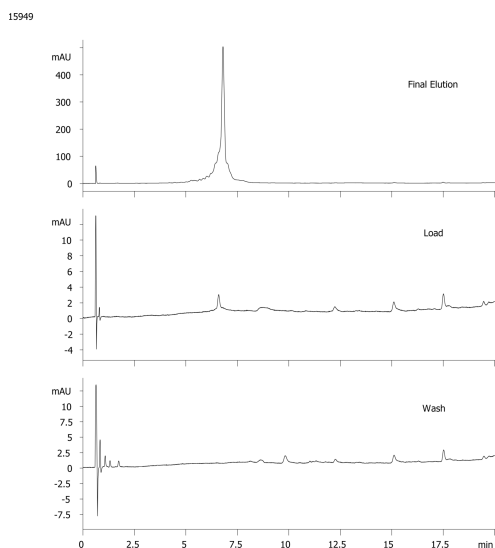

### ANALYTES:

1 Crude DNA 40nt

