

## Butter Triglycerides on ZB-5HT 15m x 0.32mm x 0.10µm

**Column:** Zebtron™ ZB-5HT, GC Cap. Column 15 m x 0.32 mm x 0.10 µm, Ea

**Phase:** 5% Phenyl 95% Dimethylpolysiloxane

**Dimensions:** 15 meters x 0.32 mm x 0.1 µm

**Order No:** 7EM-G015-02

**Oven Profile:** 100°C to 400°C at 14°C/min for 10 min

**Carrier Gas:** Constant Flow Helium, 1.8 mL/min

**Injection:** On-Column :1 2 µL @ 103°C


**Detection:** Flame Ionization (FID) (400°C)

**Analyst Note:** Recommended Accessories:  
Recommended Liner: Zebtron PLUS Liner for Agilent & Thermo  
4mm ID Direct Connect

Liner Part No.: AG2-0A50-05



Products used in this application:

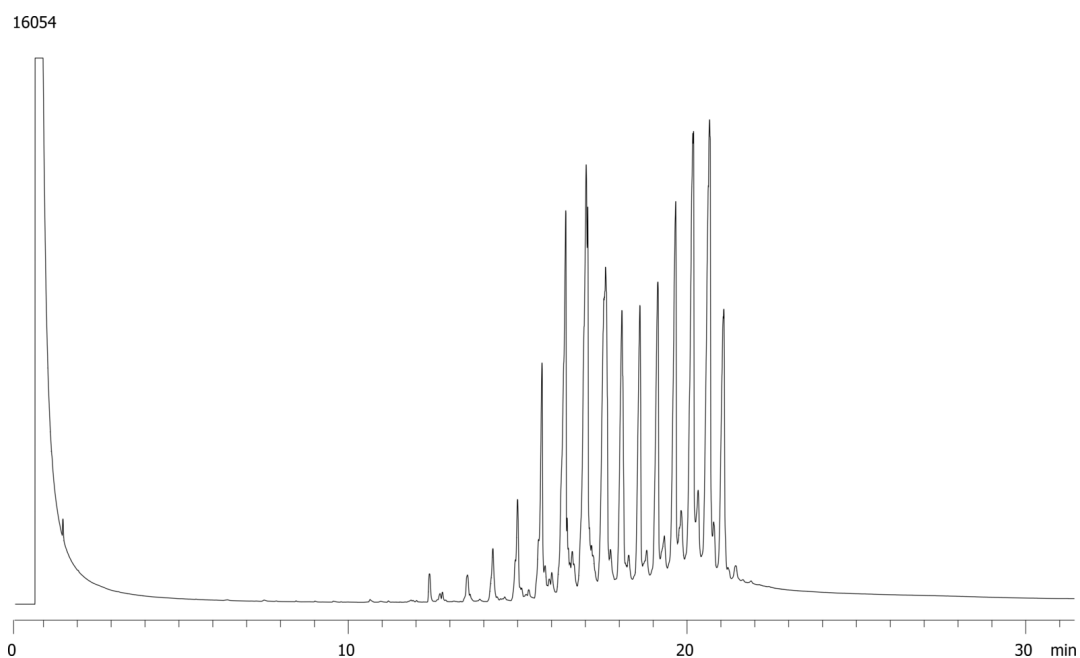

### ANALYTES:

1 butter

