

Blood Alcohols on ZB-1701

Column: Zebtron™ ZB-1701, GC Cap. Column 30 m x 0.32 mm x 0.25 µm, Ea

Phase: 14% Cyanopropylphenyl 86% Dimethylpolysiloxane

Dimensions: 30 meters x 0.32 mm x 0.25 µm

Order No: 7HM-G006-11

Oven Profile: 30°C Isothermal

Carrier Gas: Constant Flow Helium, 1.8 mL/min


Injection: Split 51:1 200 µL @ 200°C

Detection: Flame Ionization (FID) (220°C)

Analyst Note: Static Headspace: 1mL 50ppm sample in 20mL vial at 70°C for 15min.



Products used in this application:

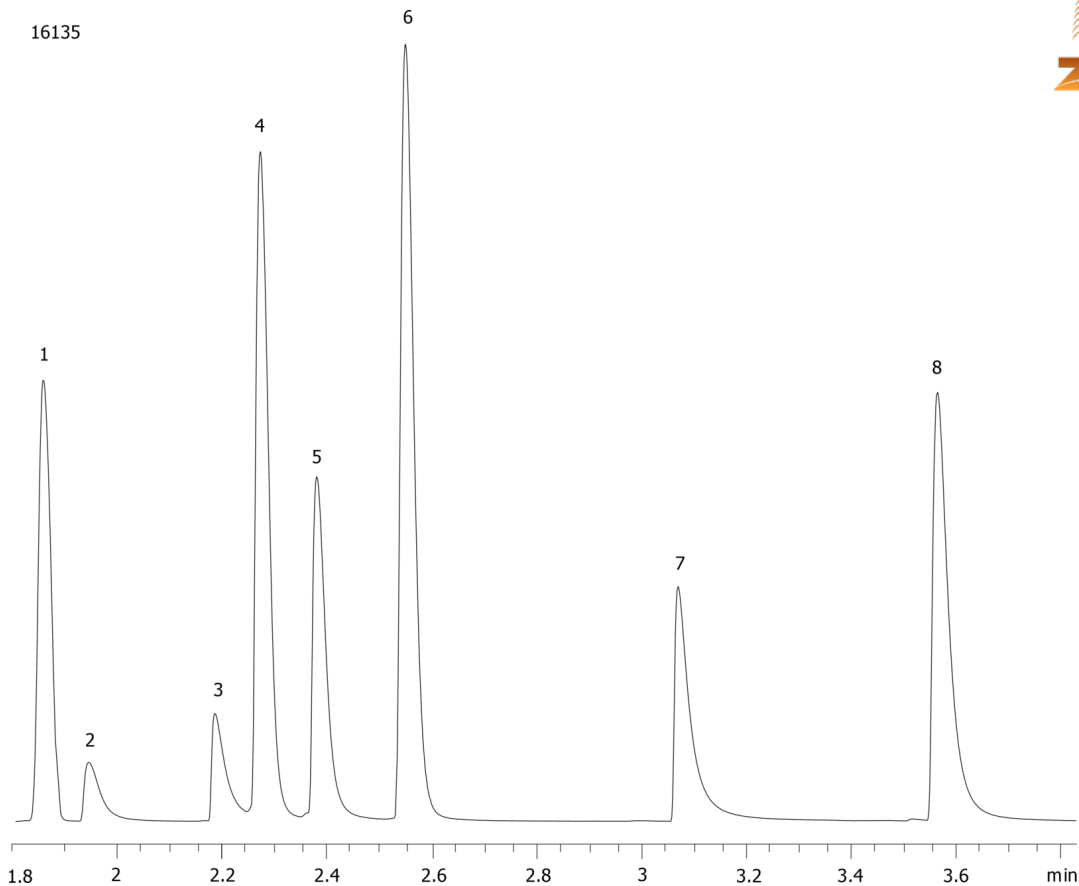

ANALYTES:

- 1 Acetaldehyde
- 2 Methanol
- 3 Ethanol
- 4 Acetone
- 5 Isopropanol
- 6 t-Butanol
- 7 n-Propanol
- 8 2-Butanol

