
## GC Application ID No.: 16165



## Low Breakdown of DDT on ZB-MR1

Zebron<sup>™</sup> ZB-MultiResidue<sup>™</sup>-1, GC Cap. Column 30 m x 0.32 mm x 0.50 µm, Ea Column: Phase: **Dimensions:** 30 meters x 0.32 mm x 0.5 µm **Order No:** 7HM-G016-17 **Oven Profile:** 100 °C for 0.5 min to 340 °C at 15 °C/min for 2 min **Carrier Gas:** Constant Flow Helium, 1.6 mL/min Injection: Splitless :1 1 µL @ 250°C **Detection:** Electron Capture (ECD) (350°C) Products used in this application: Analyst Note: Analyte was 250 pg on column. 16165

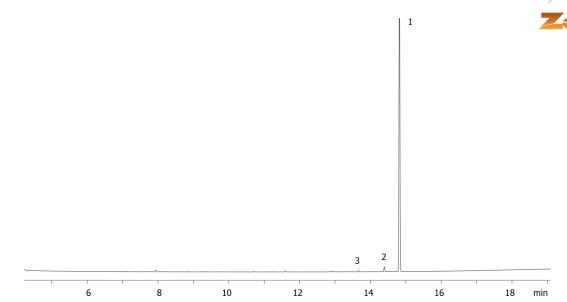

## **ANALYTES:**

- DDT 1
- DDD 2
- 3 DDE

©2025 Phenomenex Inc. All rights reserved.

Phenomenex products are available worldwide.

www.phenomenex.com

For more information contact your Phenomenex Representative at support@phenomenex.com

support@phenomenex.com