

## Nitrogen and Phosphorous Pesticides on ZB-MR1

**Column:** Zebtron™ ZB-MultiResidue™-1, GC Cap. Column 30 m x 0.25 mm x 0.25 μm, Ea

**Phase:** Unknown or Proprietary

**Dimensions:** 30 meters x 0.25 mm x 0.25 μm

**Order No:** 7HG-G016-11

**Oven Profile:** 60 °C for 0.5 min to 320 °C @ 10 °C/min for 5 min

**Carrier Gas:** Constant Flow Helium, 0.9 mL/min


**Injection:** Splitless :1 1 μL @ 260°C

**Detection:** Mass Selective (MSD) (180°C)

16377



Products used in this application:

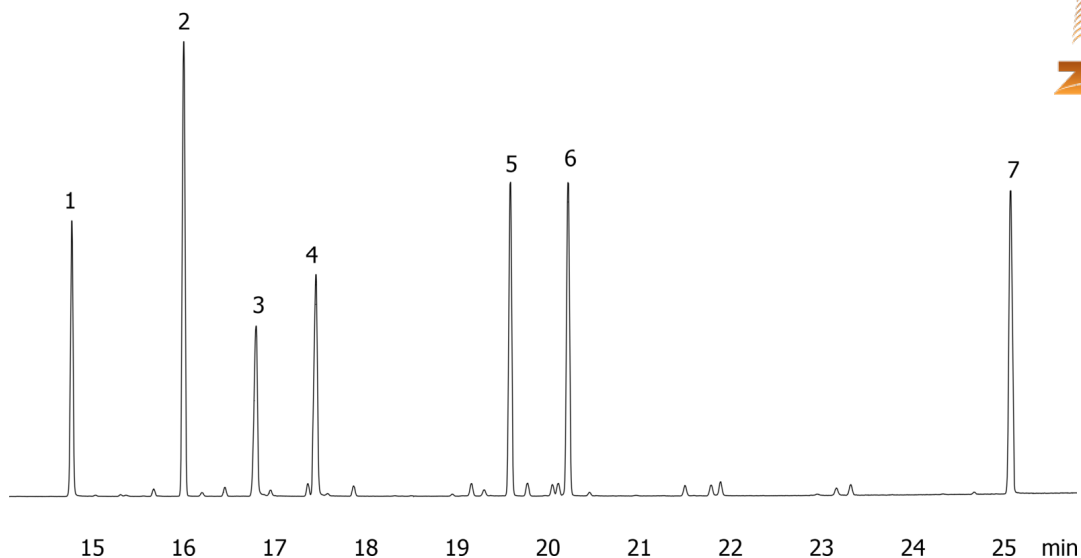

### ANALYTES:

- 1 TEPP
- 2 Sulfotep
- 3 Monochrotophos
- 4 Dimethoate
- 5 Malathion
- 6 Parathion
- 7 EPN

