

# HPLC Application

ID No.: 16398

## Oligo WAX - AP-A&G sites & DNA 20mer

**Column:** N/A  
**Dimensions:** N/A  
**Order No:** N/A  
**Elution Type:** N/A  
**Eluent A:** Water  
**Eluent B:** Acetonitrile  
**Eluent C:** 100mM Tris pH 8  
**Eluent D:** 2M NaCl

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B	Pct C	Pct D
	1	0	70	10	20	0
	2	20	10	10	20	60

**Flow Rate:** 2.2 mL/min


**Col. Temp.:** ambient

**Detection:** UV-Vis Abs.-Diode Array (PDA) @ (ambient)

**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Clarity Oligo-WAX 4 x 3.0mm ID, 10/Pk Part No.: AJ0-8324

- Holder Part No.: KJ0-4282

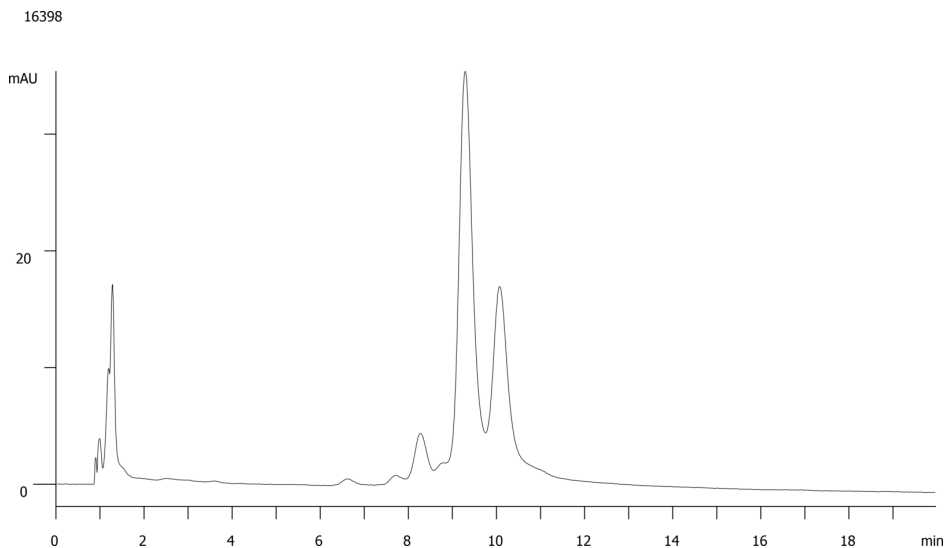

### ANALYTES:

- 1 Depurinated A&G & 20mer DNA

