


## Organophosphorus Pesticides on ZB-XLB

**Column:** Zebtron™ ZB-XLB, GC Cap. Column 30 m x 0.25 mm x 0.25 µm, Ea  
**Phase:** Proprietary XLB Phase  
**Dimensions:** 30 meters x 0.25 mm x 0.25 µm  
**Order No:** 7HG-G019-11  
**Oven Profile:** 100 °C to 200 °C @ 20 °C/min to 220 °C @ 5°C/min to 300 °C @ 15 °C/min for 2 min.  
**Carrier Gas:** Constant Flow Helium, 1 mL/min  
**Injection:** Split 5:1 1 µL @ 250°C  
**Detection:** Mass Selective (MSD) (180°C)  
**Analyst Note:** Analytes 25ppm in Isooctane



Products used in this application:

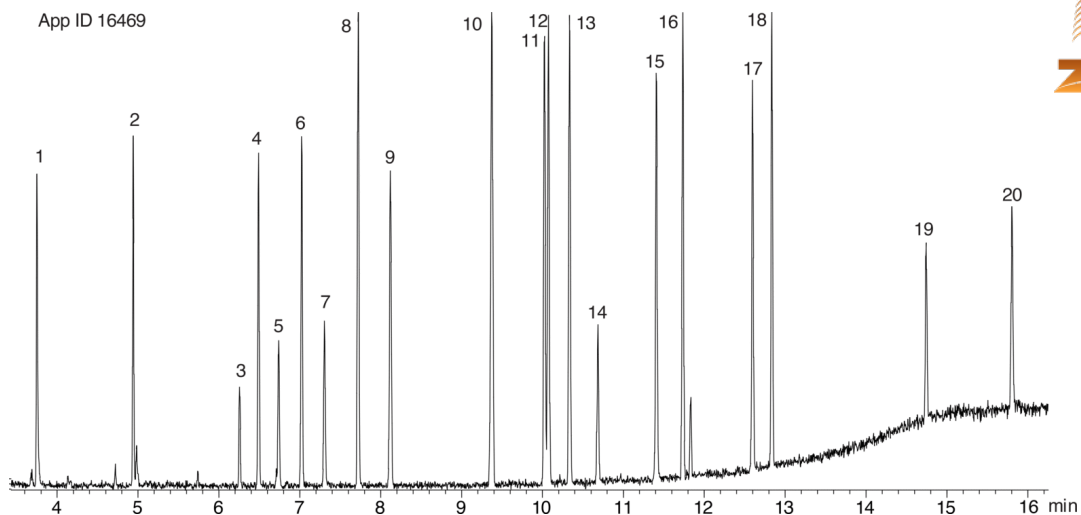

## Organophosphorus Pesticides on ZB-XLB

### ANALYTES:

- 1 Dichlorvos
- 2 Mevinphos
- 3 Demeton
- 4 Ethoprop
- 5 Naled
- 6 Phorate
- 7 Demeton
- 8 Diazinon
- 9 Disulfoton
- 10 Ronnel
- 11 Chlorpyrifos
- 12 Fenthion
- 13 Trichloronate
- 14 Merphos
- 15 Tetrachlorvinphos
- 16 Tokuthion
- 17 Fensulfothion
- 18 Bolstar
- 19 Azinphos-methyl
- 20 Coumaphos

