

## Glycols and Benzene in Toothpaste Standard on ZB-WAXplus

**Column:** Zebron™ ZB-WAXplus™, GC Cap. Column 30 m x 0.25 mm x 0.25 µm, Ea

**Phase:** Polyethylene Glycol (PEG)

**Dimensions:** 30 meters x 0.25 mm x 0.25 µm

**Order No:** 7HG-G013-11

**Oven Profile:** 80 °C to 180 °C @ 20 °C/min for 20 min.

**Carrier Gas:** Constant Flow Helium, 3 mL/min


**Injection:** Split 20:1 0.5 µL @ 225°C

**Detection:** Flame Ionization (FID) (250°C)

**Analyst Note:** Analytes at 0.5 mg/ml



Products used in this application:

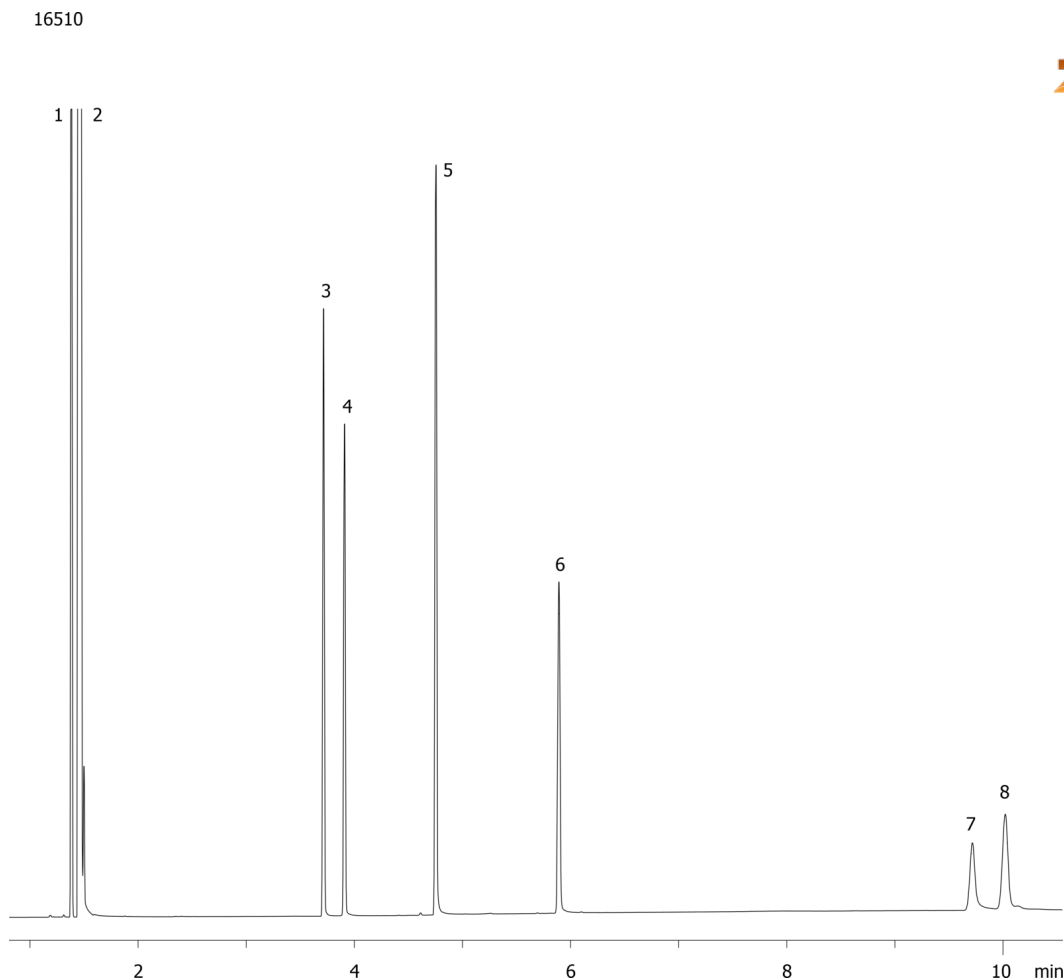

### ANALYTES:

- 1 Benzene
- 2 Acetonitrile
- 3 1,2-Propanediol
- 4 Ethylene glycol
- 5 1,3-Propanediol
- 6 Diethylene glycol
- 7 Glycerol
- 8 Triethylene glycol

© 2014 Phenomenex, Inc. All rights reserved.

For more information contact your Phenomenex Representative at [support@phenomenex.com](mailto:support@phenomenex.com)



Phenomenex products are available worldwide.

[www.phenomenex.com](http://www.phenomenex.com)

[support@phenomenex.com](mailto:support@phenomenex.com)