

## Halogenated Pesticides Mix # 2 on ZB-MR1

**Column:** Zebtron™ ZB-MultiResidue™-1, GC Cap. Column 30 m x 0.25 mm x 0.25 μm, Ea

**Phase:** Proprietary Pesticides Phase

**Dimensions:** 30 meters x 0.25 mm x 0.25 μm

**Order No:** 7HG-G016-11

**Oven Profile:** 40 °C for 1 min to 100 °C @ 15 °C/min to 225 °C @ 5 °C/min to 320 °C @ 15 °C/min for 5 min

**Carrier Gas:** Constant Flow Helium, 0.7 mL/min


**Injection:** On-Column :1 0.5 μL @ 43°C

**Detection:** Mass Selective (MSD) (335°C)

**Analyst Note:** Analytes 10 ppm in Acetone



Products used in this application:

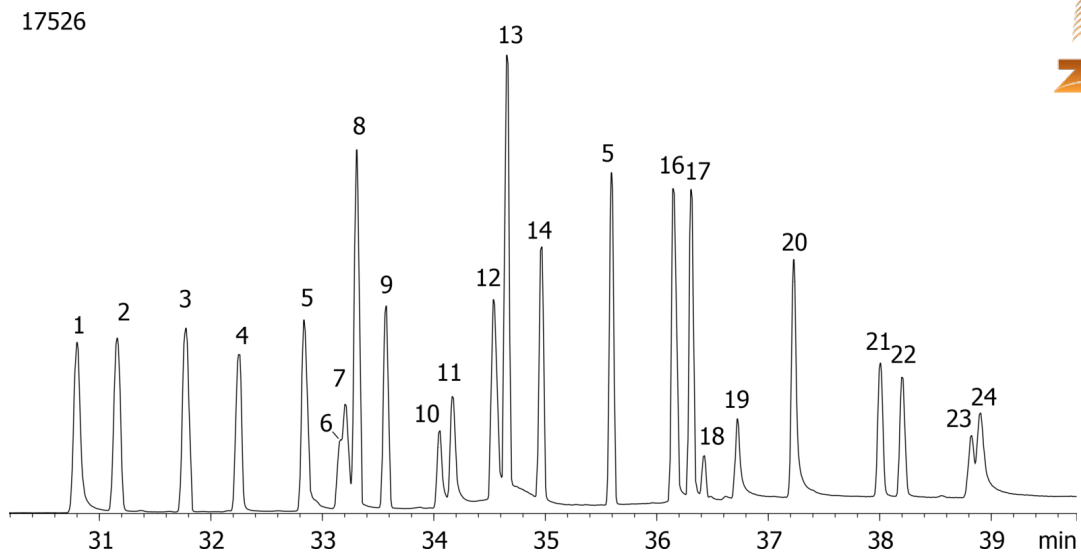

## Halogenated Pesticides Mix # 2 on ZB-MR1

### ANALYTES:

- 1 Penconazole
- 2 Procymidone
- 3 Endosulfan I
- 4 Fludioxonil
- 5 Myclobutanil
- 6 Cyproconazole
- 7 Cyproconazole
- 8 o,p'-DDT
- 9 Endosulfan II
- 10 Propiconazole
- 11 Propiconazole
- 12 Tebuconazole
- 13 Bifenthrin
- 14 Bromopropylate
- 15 Tetradifon
- 16 Fenarimol
- 17 cis-Permethrin
- 18 trans-Permethrin
- 19 Prochloraz
- 20 Fenbuconazole
- 21 Fenvalerate
- 22 Fenvalerate
- 23 Difenoconazole
- 24 Difenoconazole

