

Polynuclear Aromatic Hydrocarbons (PAHs) on ZB-MR1

Column: Zebron™ ZB-MultiResidue™-1, GC Cap. Column 30 m x 0.25 mm x 0.25 μm, Ea

Phase: Proprietary Pesticides Phase

Dimensions: 30 meters x 0.25 mm x 0.25 μm

Order No: 7HG-G016-11

Oven Profile: 40 °C for 1 min to 100 °C @ 15 °C/min to 225 °C @ 5 °C/min to 320 °C @ 15 °C/min for 5 min

Carrier Gas: Constant Flow Helium, 0.7 mL/min

Injection: On-Column :1 0.5 μL @ 43°C

Detection: Mass Selective (MSD) (335°C)


Analyst Note: Analytes 10 ppm in Dichloromethane



Products used in this application:



17557

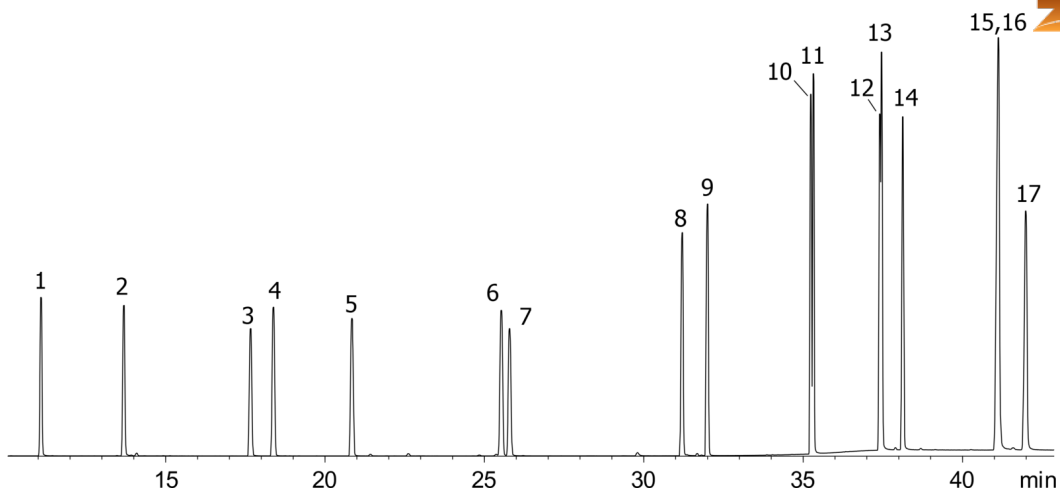

Polynuclear Aromatic Hydrocarbons (PAHs) on ZB-MR1

ANALYTES:

- 1 Naphthalene
- 2 2-Methylnaphthalene
- 3 Acenaphthylene
- 4 Acenaphthene
- 5 Fluorene
- 6 Phenanthrene
- 7 Anthracene
- 8 Fluoranthene
- 9 Pyrene
- 10 Benz[a]anthracene
- 11 Chrysene
- 12 Benzo[b]fluoranthene
- 13 Benzo[k]fluoranthene
- 14 Benzo[a]pyrene
- 15 Dibenz[a,h]anthracene
- 16 Indeno[1,2,3-cd]pyrene
- 17 Benzo[g,h,i]perylene

