

Nitrogen and Phosphorous Pesticides (EPA 507) Mix # 6 on ZB-MR2

Column: Zebtron™ ZB-MultiResidue™-2, GC Cap. Column 30 m x 0.25 mm x 0.20 µm, Ea

Phase: Proprietary Pesticides Phase

Dimensions: 30 meters x 0.25 mm x 0.2 µm

Order No: 7HG-G017-10

Oven Profile: 40 °C for 1 min to 100 °C @ 15 °C/min to 225 °C @ 5 °C/min to 320 °C @ 15 °C/min for 5 min

Carrier Gas: Constant Flow Helium, 0.8 mL/min

Injection: On-Column :1 0.5 µL @ 43°C

Detection: Mass Selective (MSD) (335°C)


Analyst Note: Analytes 10 ppm in Acetone



Products used in this application:



17593

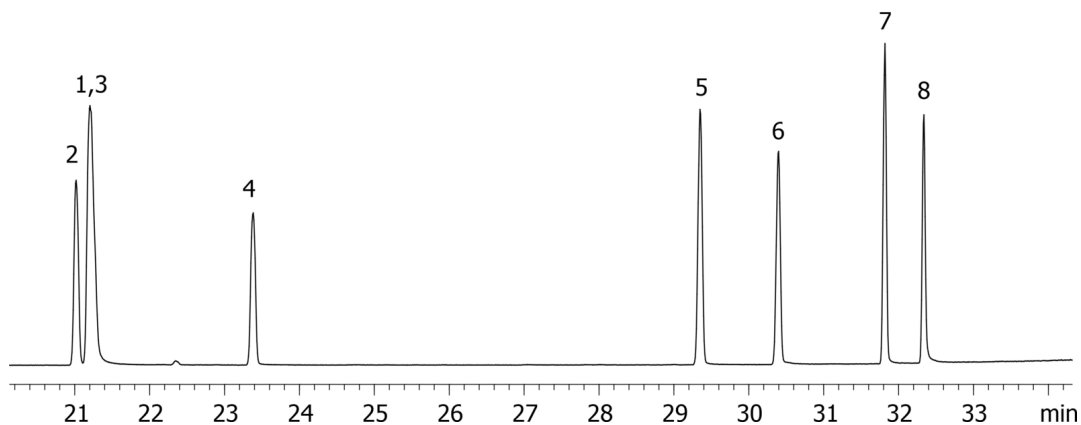

ANALYTES:

- 1 Propachlor
- 2 Trifluralin
- 3 Benefin
- 4 Profluralin
- 5 Isopropalin
- 6 Pendimethalin
- 7 Oxadiazon
- 8 Oxyfluorfen

