

## EPA 8270 Internal Standard Mix on ZB-MR2

**Column:** Zebtron<sup>™</sup> ZB-MultiResidue<sup>™</sup>-2, GC Cap. Column 30 m x 0.25 mm x 0.20 μm, Ea

**Phase:** Proprietary Pesticides Phase

**Dimensions:** 30 meters x 0.25 mm x 0.2 μm

**Order No:** 7HG-G017-10

**Oven Profile:** 40 °C for 1 min to 100 °C @ 15 °C/min to 225 °C @ 5 °C/min to 320 °C @ 15 °C/min for 5 min

**Carrier Gas:** Constant Flow Helium, 0.8 mL/min

**Injection:** On-Column :1 0.5 μL @ 43°C

**Detection:** Mass Selective (MSD) (335°C)


**Analyst Note:** Analytes 10 ppm in Acetone



Products used in this application:



17596

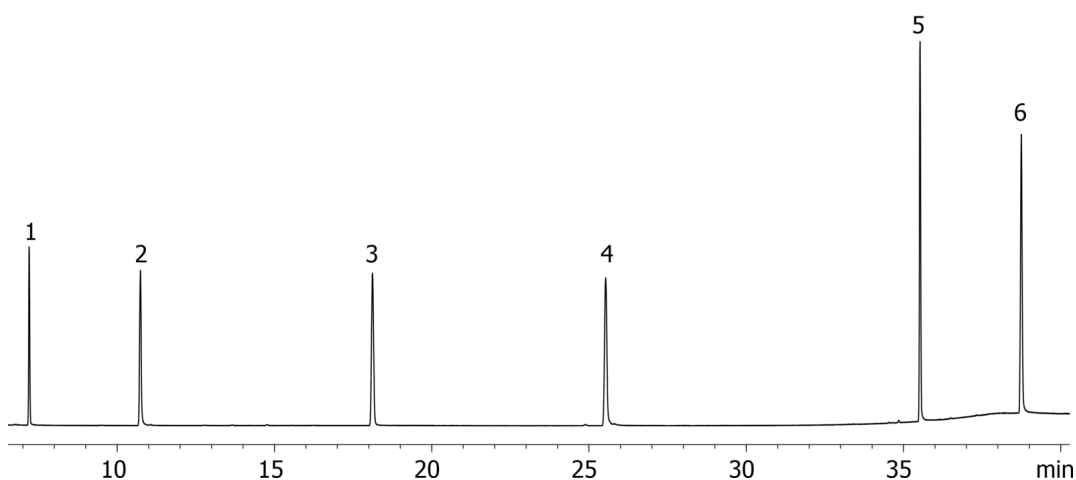

### ANALYTES:

- 1 1,4-Dichlorobenzene-d4
- 2 Naphthalene-d8
- 3 Acenaphthene-d10
- 4 Phenanthrene-d10
- 5 Chrysene-d12
- 6 Perylene-d12

