

# HPLC Application

ID No.: 17691

## Axia Packed Gemini 5 µm C18 (low pH loading - Diphenhydramine and Propranolol - 104mg/40µL)

**Column:** Gemini® 5 µm C18 110 Å, LC Column 100 x 21.2 mm, AXIA™ Packed, Ea

**Dimensions:** 100 x 21.2 mm ID

**Order No:** 00D-4435-P0-AX

**Elution Type:** Gradient

**Eluent A:** 0.5 % TFA in Water

**Eluent B:** Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	95	5
	2	5	5	95

**Flow Rate:** 30 mL/min

**Col. Temp.:** ambient

**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)

**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard PREP Cartridge, Gemini C18 15 x 21.2mm ID, Ea Part No.: AJ0-7846


- Holder Part No.: AJ0-8223



Products used in this application:



17691

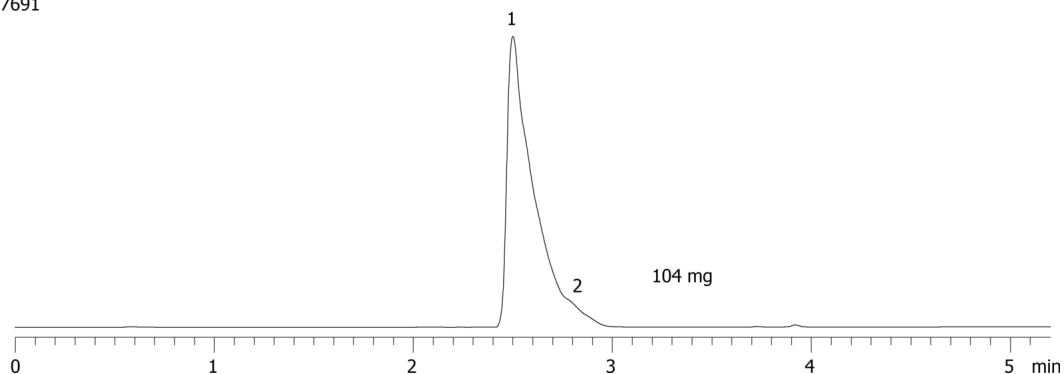

### ANALYTES:

1 Diphenhydramine

2 Propranolol

