

# HPLC Application

ID No.: 17721

## Axia Packed Gemini 5 µm C18 (high pH loading - 3 component mix - 106mg/140µL)

**Column:** Gemini<sup>®</sup> 5 µm C18 110 Å, LC Column 50 x 21.2 mm, AXIA<sup>™</sup> Packed, Ea

**Dimensions:** 50 x 21.2 mm ID

**Order No:** 00B-4435-P0-AX

**Elution Type:** Gradient

**Eluent A:** 0.2 % NH<sub>4</sub>OH in Water

**Eluent B:** Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	95	5
	2	5	5	95

**Flow Rate:** 30 mL/min

**Col. Temp.:** ambient

**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)

**Analyst Note:** SecurityGuard<sup>™</sup> Guard Cartridge System extends column lifetime.

- SecurityGuard PREP Cartridge, Gemini C18 15 x 21.2mm ID, Ea Part No.: AJ0-7846


- Holder Part No.: AJ0-8223



Products used in this application:



17721

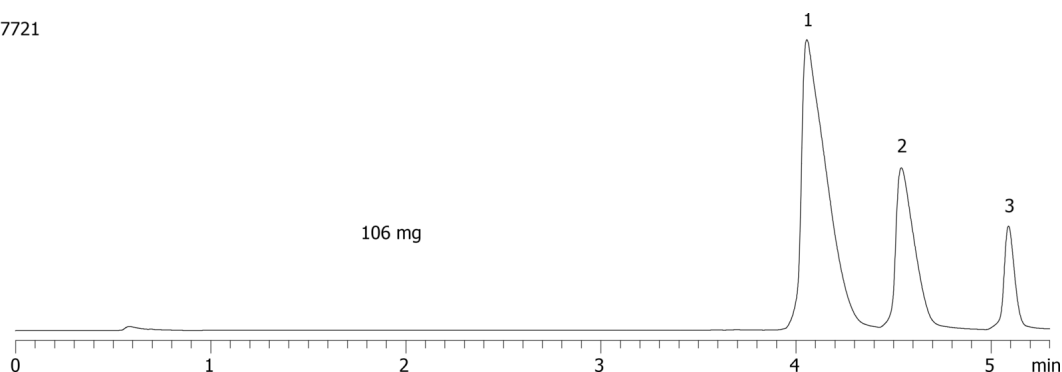

### ANALYTES:

- 1 Diphenhydramine
- 2 Oxybutynin
- 3 Terfenadine

