

## Organophosphorous Pesticides Mix 4 on ZB-50


**Column:** Zebron<sup>™</sup> ZB-50, GC Cap. Column 30 m x 0.25 mm x 0.25  $\mu$ m, Ea  
**Phase:** 50% Phenyl 50% Dimethylpolysiloxane  
**Dimensions:** 30 meters x 0.25 mm x 0.25  $\mu$ m  
**Order No:** 7HG-G004-11  
**Oven Profile:** 40 °C for 1 min to 100 °C @ 15 °C/min to 225 °C @ 5 °C/min to 320 °C @ 15 °C/min for 8 min  
**Carrier Gas:** Constant Flow Helium, 0.7 mL/min  
**Injection:** On-Column :1 0.5  $\mu$ L @ 43°C  
**Detection:** Mass Selective (MSD) (335°C)  
**Analyst Note:** Analytes are 10 ppm in Acetone



Products used in this application:



17815

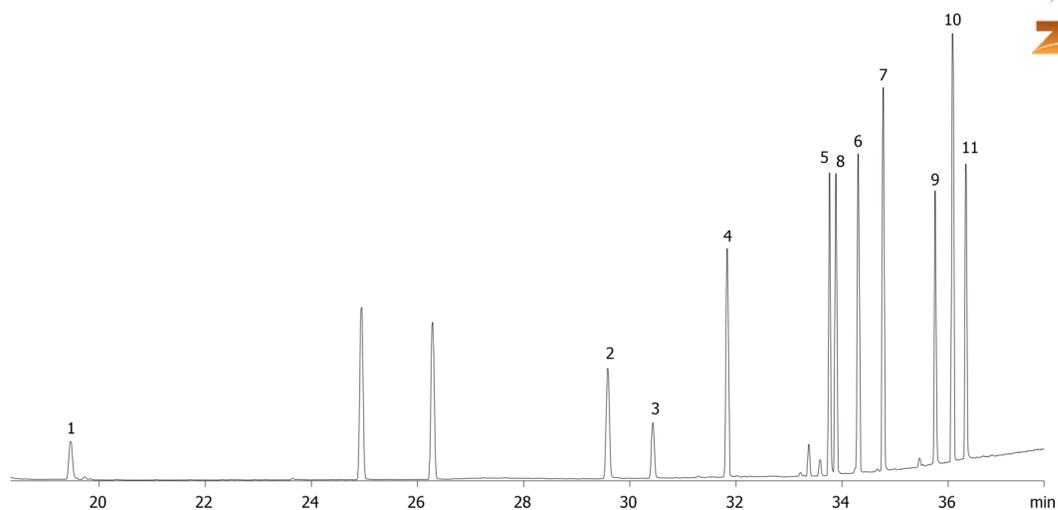

## Organophosphorous Pesticides Mix 4 on ZB-50

### ANALYTES:

- 1 Trichlorfon
- 2 Parathion OA
- 3 Chlorpyrifos OA
- 4 Mecarbam
- 5 Chlorthiophos
- 6 Carbophenothion
- 7 Cyanofenphos
- 8 Tris(2-butoxyethyl) phosphate
- 9 Pyridaphenthion
- 10 Leptophos
- 11 Pyrazophos

