

# GC Application


ID No.: 17840

## Nitrogen and Phosphorous Pesticides (EPA 507) Mix # 3 on ZB-50

**Column:** Zebtron™ ZB-50, GC Cap. Column 30 m x 0.25 mm x 0.25 µm, Ea  
**Phase:** 50% Phenyl 50% Dimethylpolysiloxane  
**Dimensions:** 30 meters x 0.25 mm x 0.25 µm  
**Order No:** 7HG-G004-11  
**Oven Profile:** 40 °C for 1 min to 100 °C @ 15 °C/min to 225 °C @ 5 °C/min to 320 °C @ 15 °C/min for 8 min  
**Carrier Gas:** Constant Flow Helium, 0.7 mL/min  
**Injection:** On-Column :1 0.5 µL @ 43°C  
**Detection:** Mass Selective (MSD) (335°C)  
**Analyst Note:** Analytes are 10 ppm in Acetone



Products used in this application:

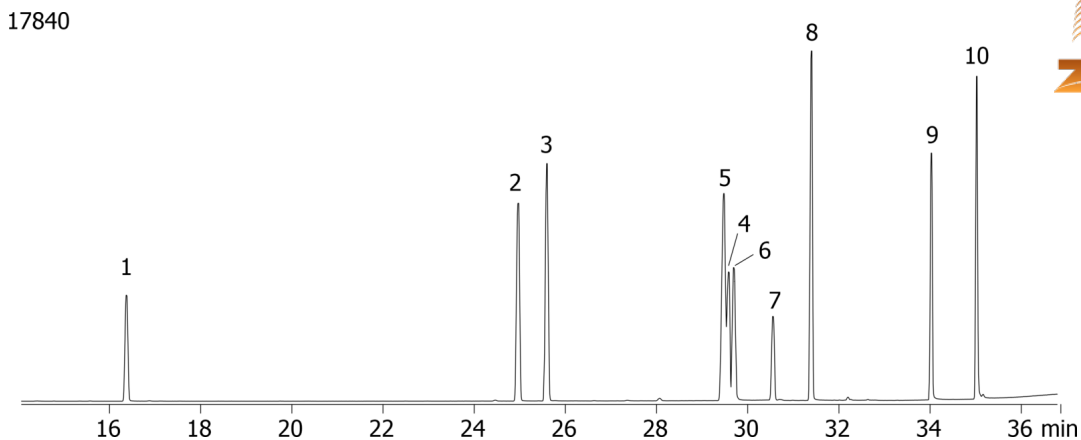

## Nitrogen and Phosphorous Pesticides (EPA 507) Mix # 3 on ZB-50

### ANALYTES:

- 1 Vernolate
- 2 Terbufos
- 3 Diazinon
- 4 Metribuzin
- 5 Metolachlor
- 6 MGK-264
- 7 MGK-264
- 8 Butachlor
- 9 Carboxin
- 10 Norflurazon

