

Benzodiazepines by MSD on a ZB-5MSi

Column: Zebtron™ ZB-5MSi, GC Cap. Column 10 m x 0.18 mm x 0.18 μm, Ea

Phase: 5% Phenyl 95% Dimethylpolysiloxane

Dimensions: 10 meters x 0.18 mm x 0.18 μm

Order No: 7CD-G018-08

Oven Profile: 200 °C to 340 °C @ 25 °C/min for 1 min

Carrier Gas: Constant Flow Helium, 1 mL/min

Injection: Split 10:1 1 μL @ 260°C

Detection: Mass Selective (MSD) (320°C)


Analyst Note: Analytes are 24ppm derivatized with MTBSTFA

Sim ions:

(Quant ion*) D5-Nordiazepam (*332,334): Nordiazepam (*327,328,383): D5-Oxazepam (*462,519): Oxazepam (*457,513,459): D5-Temazepam (*362,288,390): Temazepam (*357,359,385): Lorazepam (*491,513,493):



Products used in this application:

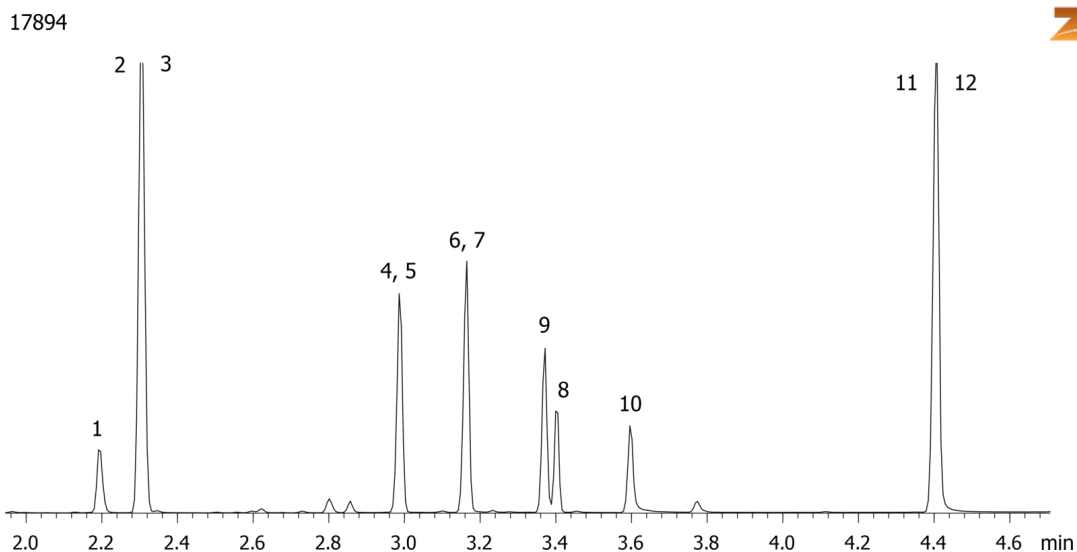

Benzodiazepines by MSD on a ZB-5MSi

ANALYTES:

- 1 Diazepam
- 2 Nordiazepam-d5
- 3 Nordiazepam
- 4 Oxazepam-d5
- 5 Oxazepam
- 6 Temazepam-d5
- 7 Temazepam
- 8 Clonazepam
- 9 Lorazepam
- 10 Alprazolam
- 11 alpha-Hydroxyalprazolam-d5
- 12 alpha-Hydroxyalprazolam

