

## Blood Alcohols on ZB-BAC-2

**Column:** Zebtron<sup>™</sup> ZB-BAC-2, GC Cap. Column 30 m x 0.32 mm x 1.20 μm, Ea

**Phase:** Proprietary Blood Alcohol Phase

**Dimensions:** 30 meters x 0.32 mm x 1.2 μm

**Order No:** 7HM-G022-25

**Oven Profile:** 40 °C Isothermal

**Carrier Gas:** Constant Flow Helium, 80 cm/sec

**Injection:** Split 5:1 0.5 μL @ 140°C


**Detection:** Flame Ionization (FID) (220°C)

**Analyst Note:** 0.050% Blood Alcohols Components in Water (1mL in 20mL Headspace Vial)

17898



Products used in this application:

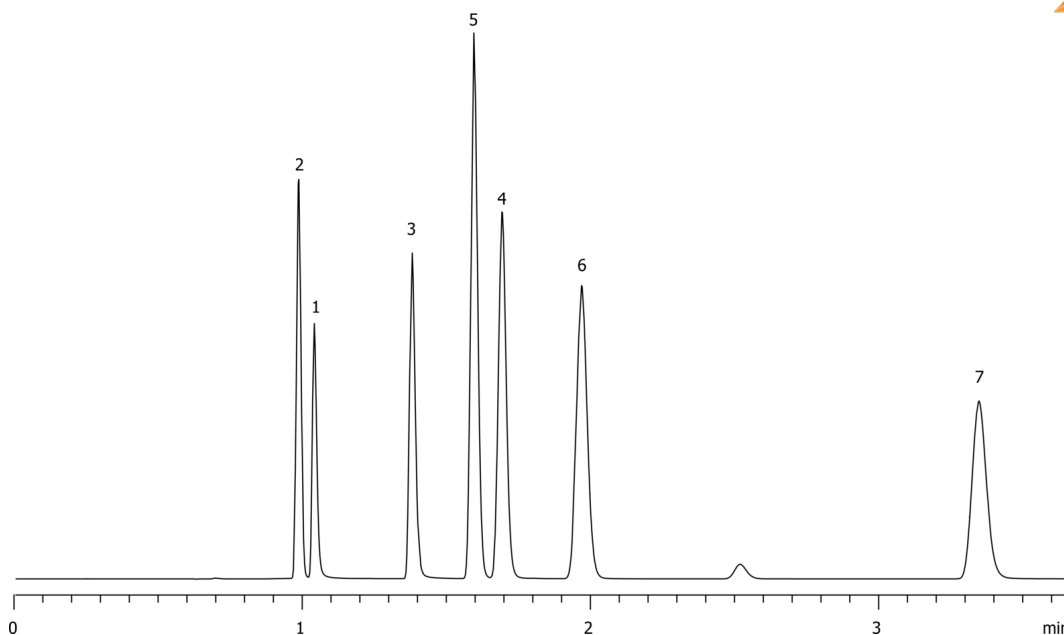

### ANALYTES:

- 1 Methanol
- 2 Acetaldehyde
- 3 Ethanol
- 4 Isopropanol
- 5 Acetone
- 6 t-Butanol
- 7 2-Butanol

