

# HPLC Application

ID No.: 17913

## EPA Method 8330: Explosive Mix Using Luna Phenyl Hexyl

**Column:** Luna® 5 µm Phenyl-Hexyl 100 Å, LC Column 250 x 4.6 mm, Ea

**Dimensions:** 250 x 4.6 mm ID

**Order No:** 00G-4257-E0

**Elution Type:** Gradient

**Eluent A:** Methanol

**Eluent B:** Water

<b>Gradient Profile:</b>	<b>Step No.</b>	<b>Time (min)</b>	<b>Pct A</b>	<b>Pct B</b>
	<b>1</b>	0	55	45
	<b>2</b>	2	55	45
	<b>3</b>	22	75	25
	<b>4</b>	22.01	55	45
	<b>5</b>	25	55	45

**Flow Rate:** 1 mL/min

**Col. Temp.:** 38 °C

**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 220 nm (ambient)

**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Phenyl 4 x 3.0mm, 10/Pk Part No.: AJ0-4351


- Holder Part No.: KJ0-4282



Products used in this application:



17913

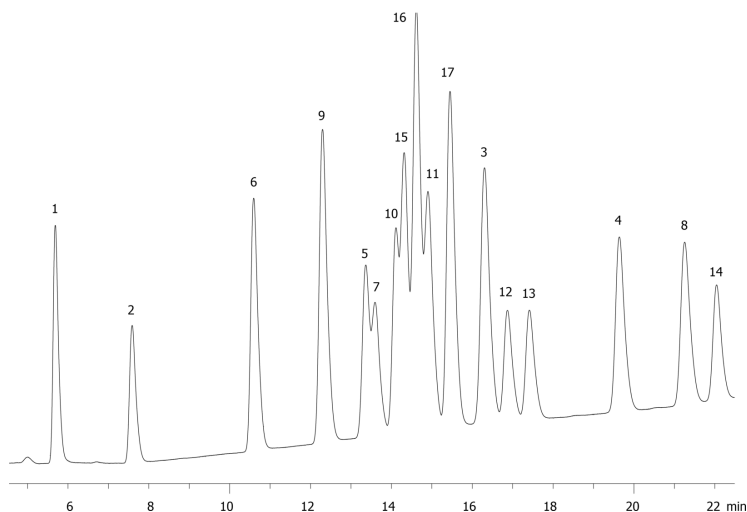

## EPA Method 8330: Explosive Mix Using Luna Phenyl Hexyl

### ANALYTES:

- 1** Octahydro-1,3,5,7-tetranitro-1,3,5,7-tetrazocine (HMX)
- 2** Hexahydro-1,3,5-trinitro-1,3,5-triazine
- 3** 1,3,5-Trinitrobenzene
- 4** Methyl-2,4,6-trinitrophenylnitramine (Tetryl)
- 5** 1,3-Dinitrobenzene
- 6** Nitrobenzene
- 7** Nitroglycerin
- 8** 2,4,6-Trinitrotoluene
- 9** 3,5-Dinitroaniline
- 10** 4-Amino-2,6 dinitrotoluene
- 11** 2 Amino-4,6 dinitrotoluene
- 12** 2,6-Dinitrotoluene
- 13** 2,4-Dinitrotoluene
- 14** Pentaerythritol tetranitrate
- 15** 2-Nitrotoluene
- 16** 4-Nitrotoluene
- 17** 3-Nitrotoluene

