

ABN Mix on Synergi 4u Max-RP 150x4.6

Column: Synergi™ 4 µm Max-RP 80 Å, LC Column 150 x 4.6 mm, Ea
Dimensions: 150 x 4.6 mm ID
Order No: 00F-4337-E0
Elution Type: Gradient
Eluent A: Phosphate Buffer, pH 2.5
Eluent B: Acetonitrile
Eluent C: 0
Eluent D: 0



Gradient Profile:	Step No.	Time (min)	Pct A	Pct B	Pct C	Pct D
	1	0	95	5	0	0
	2	8	20	80	0	0
	3	10	20	80	0	0
	4	10.5	95	5	0	0
	5	12	95	5	0	0

Flow Rate: 1 mL/min


Col. Temp.: ambient

Detection: UV-Vis Abs.-Diode Array (PDA) @ 254 210 (ambient)

Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, C12 4 x 3.0mm, 10/Pk Part No.: AJ0-6074

- Holder Part No.: KJ0-4282

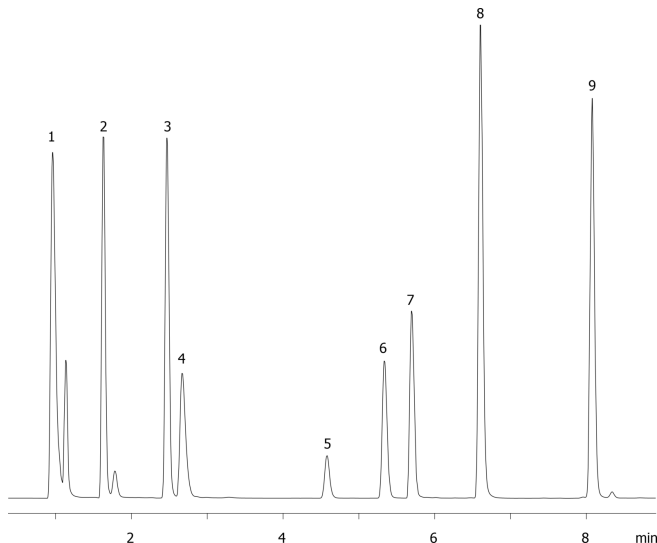

ANALYTES:

- 1** Pyridine
Retention Time: 1.1 min
- 2** Quinidine
Retention Time: 1.6 min
- 3** Sulfathiazole
Retention Time: 2.5 min
- 4** Triprolidine
Retention Time: 2.7 min
- 5** Benzyl alcohol
Retention Time: 4.6 min
- 6** Nortryptiline
Retention Time: 5.3 min
- 7** Phenol
Retention Time: 6.6 min
- 8** 3-Methyl 4-benzoic acid
Retention Time: 8.1 min
- 9** 3-Methyl salicylaldehyde
Retention Time: 8.1 min

© 2024 Phenomenex, Inc. All rights reserved.

For more information contact your Phenomenex Representative at support@phenomenex.com

Phenomenex products are available worldwide.
www.phenomenex.com

support@phenomenex.com

- 8** 3-Methyl 4-benzoic acid
Retention Time: 8.1 min
- 9** 3-Methyl salicylaldehyde
Retention Time: 8.1 min