

# HPLC Application

ID No.: 17980

## EPA Method 8330: Explosive Mix Using Synergi 2.5u Hydro Rp

**Column:** Synergi<sup>™</sup> 2.5  $\mu$ m Hydro-RP 100 Å, LC Column 100 x 2 mm, Ea

**Dimensions:** 100 x 2 mm ID

**Order No:** 00D-4387-B0

**Elution Type:** Gradient

**Eluent A:** Methanol

**Eluent B:** Water


Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	55	45
	2	0.5	55	45
	3	8	75	25
	4	10	75	25

**Flow Rate:** 0.2 mL/min

**Col. Temp.:** ambient

**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 220 nm (ambient)

17980

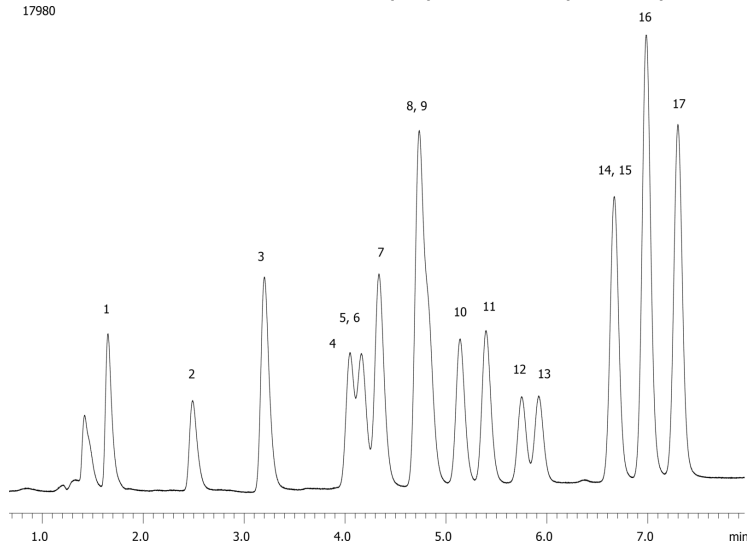

## EPA Method 8330: Explosive Mix Using Synergi 2.5u Hydro Rp

### ANALYTES:

- 1** Octahydro-1,3,5,7-tetranitro-1,3,5,7-tetrazocine (HMX)
- 2** Hexahydro-1,3,5-trinitro-1,3,5-triazine
- 3** 1,3,5-Trinitrobenzene
- 4** Methyl-2,4,6-trinitrophenylnitramine (Tetryl)
- 5** 1,3-Dinitrobenzene
- 6** Nitrobenzene
- 7** Nitroglycerin
- 8** 2,4,6-Trinitrotoluene
- 9** 3,5-Dinitroaniline
- 10** 4-Amino-2,6 dinitrotoluene
- 11** 2 Amino-4,6 dinitrotoluene
- 12** 2,6-Dinitrotoluene
- 13** 2,4-Dinitrotoluene
- 14** Pentaerythritol tetranitrate
- 15** 2-Nitrotoluene
- 16** 4-Nitrotoluene
- 17** 3-Nitrotoluene

