

## Blood Alcohols using Hydrogen on ZB-BAC-1

**Column:** Zebtron™ ZB-BAC-1, GC Cap. Column 30 m x 0.53 mm x 3.00 μm, Ea

**Phase:** Proprietary Blood Alcohol Phase

**Dimensions:** 30 meters x 0.53 mm x 3 μm

**Order No:** 7HK-G021-36

**Oven Profile:** 40 °C Isothermal

**Carrier Gas:** Constant Flow Hydrogen, 80 cm/sec

**Injection:** Split 5:1 1 μL @ 150°C

**Detection:** Flame Ionization (FID) (250°C)


**Analyst Note:** All analytes are 0.100% in water



Products used in this application:



18025

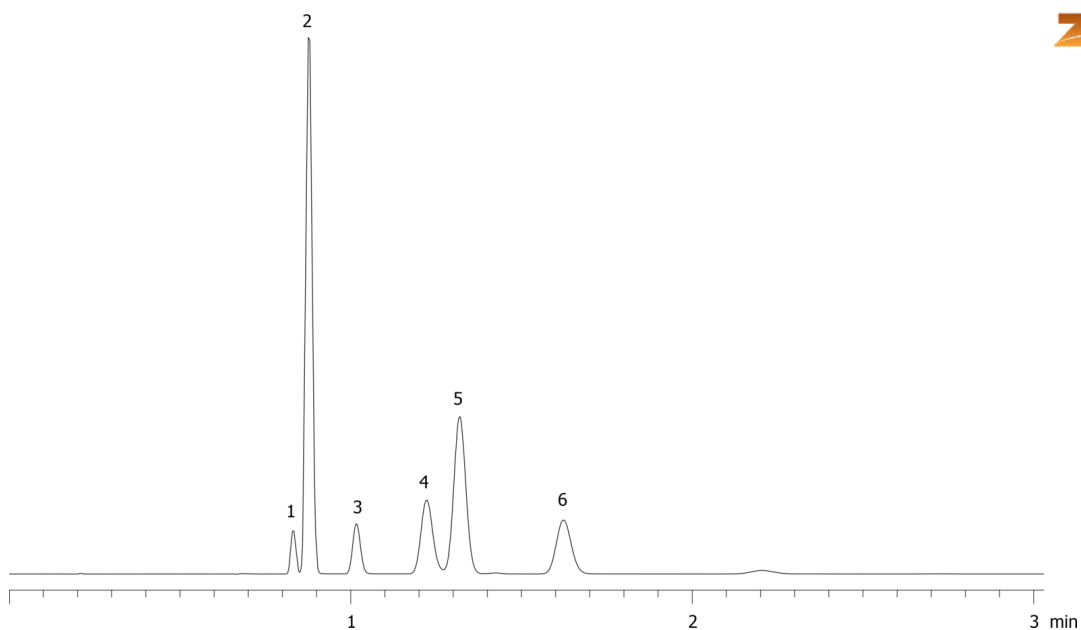

### ANALYTES:

- 1 Methanol
- 2 Acetaldehyde
- 3 Ethanol
- 4 Isopropanol
- 5 Acetone
- 6 n-Propanol

