

## Blood Alcohols on ZB-BAC-1

**Column:** Zebtron™ ZB-BAC-1, GC Cap. Column 30 m x 0.32 mm x 1.80 μm, Ea

**Phase:** Proprietary Blood Alcohol Phase

**Dimensions:** 30 meters x 0.32 mm x 1.8 μm

**Order No:** 7HM-G021-31

**Oven Profile:** 40 °C Isothermal

**Carrier Gas:** Constant Flow Helium, 80 cm/sec

**Injection:** Split 15:1 1 μL @ 150°C


**Detection:** Flame Ionization (FID) (250°C)

**Analyst Note:** Analytes and internal standards were 0.100% in water

18059



Products used in this application:

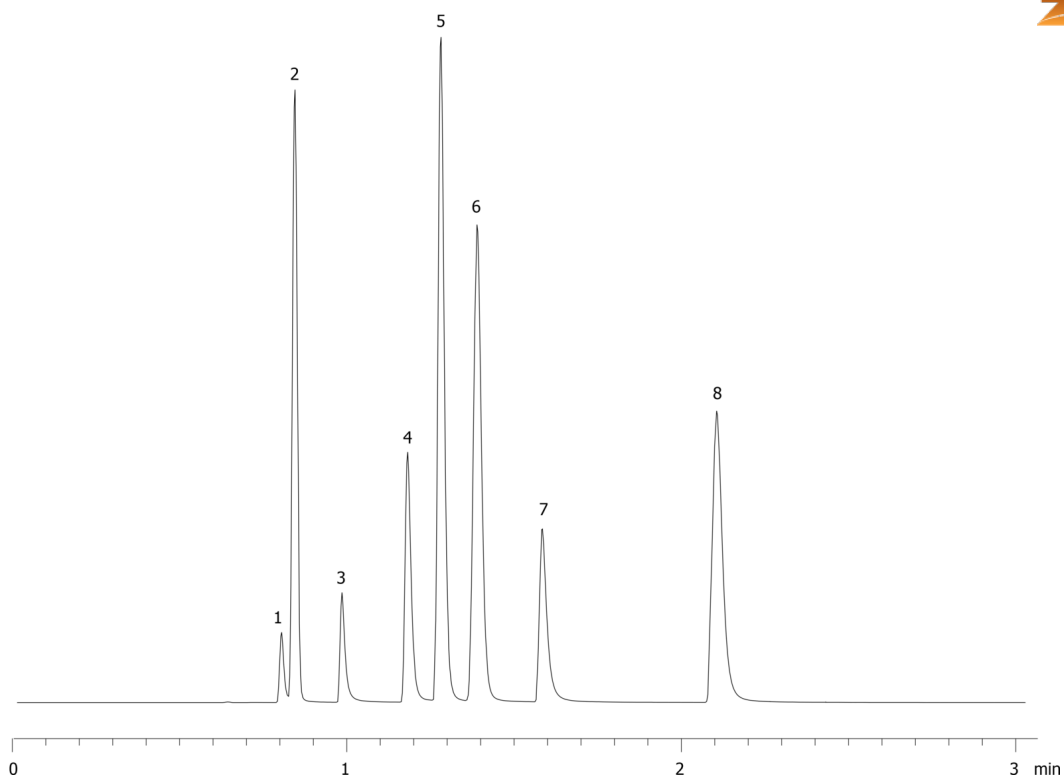

### ANALYTES:

- 1 Methanol
- 2 Acetaldehyde
- 3 Ethanol
- 4 Isopropanol
- 5 Acetone
- 6 t-Butanol
- 7 n-Propanol
- 8 2-Butanol

