

# HPLC Application

ID No.: 18213

## Trans Stilbene Oxide on a Lux 5u Cellulose 2 Axia 30mm ID column

**Column:** Lux<sup>®</sup> 5 µm Cellulose-2, LC Column 250 x 30 mm, AXIA<sup>™</sup> Packed, Ea

**Dimensions:** 250 x 30 mm ID

**Order No:** 00G-4457-U0-AX

**Elution Type:** Isocratic

**Eluent A:** Hexane

**Eluent B:** IPA

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	90	10

**Flow Rate:** 40 mL/min

**Col. Temp.:** ambient

**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 220 nm (ambient)

**Analyst Note:** SecurityGuard<sup>™</sup> Guard Cartridge System extends column lifetime.

- SecurityGuard PREP Cartridge, Lux Cellulose-2, 15 x 30mm ID, Ea Part No.: AJ0-8401


- Holder Part No.: AJ0-8277



Products used in this application:



18213

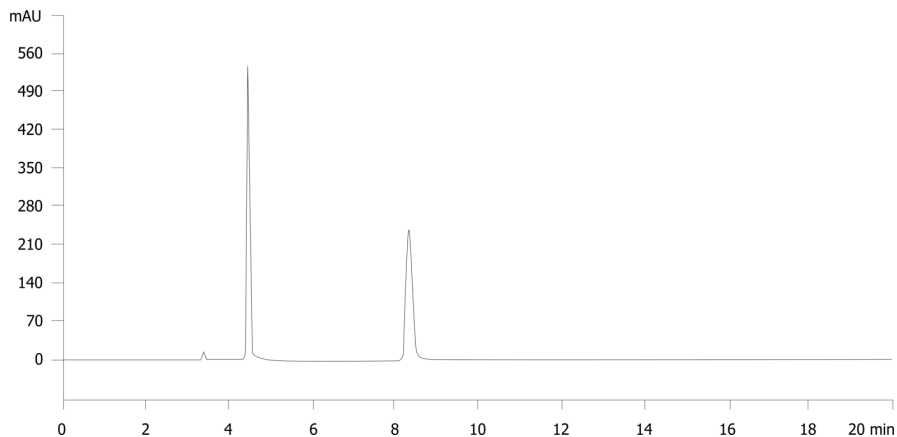

### ANALYTES:

1 trans-Stilbene oxide