


## EPA 8270 Internal Standard Mix on ZB-1701P

**Column:** Zebtron™ ZB-1701, GC Cap. Column 30 m x 0.25 mm x 0.25 μm, Ea  
**Phase:** 14% Cyanopropylphenyl 86% Dimethylpolysiloxane  
**Dimensions:** 30 meters x 0.25 mm x 0.25 μm  
**Order No:** 7HG-G006-11  
**Oven Profile:** 40 °C for 1 min to 100 °C @ 15 °C/min to 225 °C @ 5 °C/min to 290 °C @ 15 °C/min for 25 min  
**Carrier Gas:** Constant Flow Helium, 0.8 mL/min  
**Injection:** On-Column :1 0.5 μL @ 43°C  
**Detection:** Mass Selective (MSD) (335°C)  
**Analyst Note:** Analytes 10 ppm in Acetone



Products used in this application:

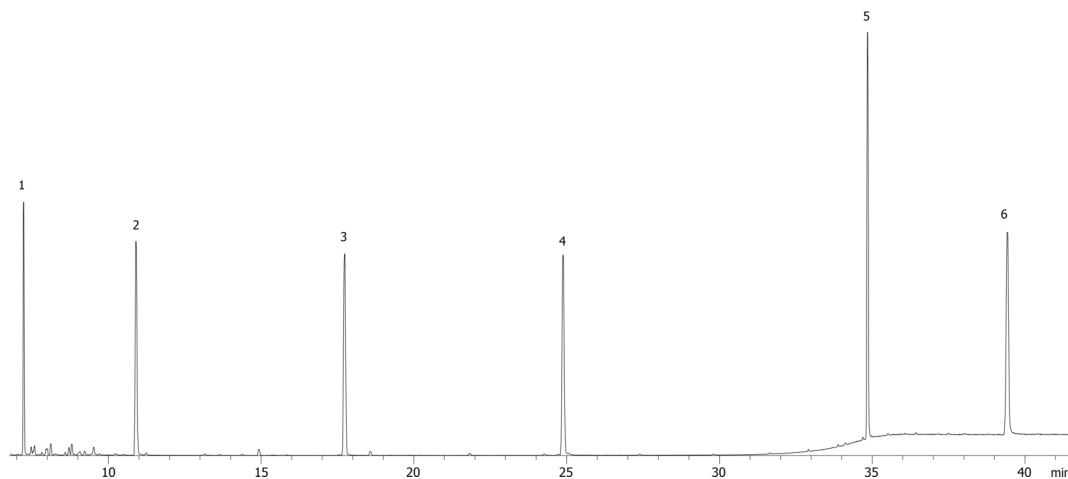

### ANALYTES:

- 1,4-Dichlorobenzene-d4
- Naphthalene-d8
- Acenaphthene-d10
- Phenanthrene-d10
- Chrysene-d12
- Perylene-d12

